L'IMP

CAANT

OCHE

EV O

VE

LEADER  MEDICA

MADE IN ITALY

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Introduction	5
The Scientific Committee	5
B.O.P.T.®	12
Ignazio Loi	12
Biological Polygon	13
Maximum Surgical Flexibility	14
B.O.P.T.® - Trust the Original	16
The Company	18
The DNA of Evolution	18
Leader Medica Worldwide	20
IL® e IL WIDE® Features	22
Two connections	25
FMC Neck	26
Straight Ahead	27
DAS	28
Decontamination	30
ATC Threads and Bodyblade®	31
Certification of the Protocol for Surface Treatments and Processes	32
Surface Morphology	34
IL® e IL WIDE®	36
IL® and IL WIDE® - Guide and coding of systems	36
Protocol and Surgical Kit	40
Surgical Protocol	40
B.O.P.T.® Vertical Technique - Kit	43
IL® and IL WIDE® Kit	44
Guided Surgery	46
Prosthetically Guided Surgery	46
Digital Workflow	47
The Software	48
Why Choose It	49
Guided Protocol and Surgical Sequence	50
Features of the IL® and IL WIDE® Guided Surgery Kit	54
Pin sleeves and guide	55
IL® and IL WIDE® Guided Surgery Kit	56
Drills Sequence Guided Surgery	58
Iconic Abutment	60
Quick and Easy Vignato	60
Iconic Box	61
Iconic Abutment	62
Iconic Abutment Prosthetic Protocol	64
Prosthetic Components	66
Coding of IL® prosthetic components	68
Unico Abutment	80
Coding of IL WIDE® prosthetic components	83
Impianto IL® Zigomatico, Nasale e Pterigoideo	92
Coding of IL® ZY/NS/PT prosthetic components	93
Surgical Recommendations	96
Bibliography	97
Certifications and Authorisations	98

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/ THE SCIENTIFIC COMMITTEE

*The **Scientific Committee** supporting Leader Medica, in its institutional and interdisciplinary composition, has played a crucial role in defining the content at the heart of this publication. This catalog therefore represents the outcome of a shared project, the result of a process of institutional and academic-scientific collaboration aimed at combining the realms of research and dissemination, in full respect of methodological rigor, intellectual integrity, and the enhancement of the cultural heritage contributed by each member of the **Scientific Committee**.*

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Brånemark's studies have effectively revolutionized dental treatment planning, thanks to proving the predictability of osseointegration, which is perfect peri-implant bone healing, predictable and maintainable over time when there is rigorous understanding of bone biology.

Undoubtedly, the future of implant research will focus on what we might call the ideal integration of supracrestal soft tissues with the implant-prosthetic unit.

All this is due to the importance of this architecture in maintaining the long-term stability of implant-prosthetic rehabilitation in order to prevent infectious and inflammatory complications and to preserve the aesthetics and the mimicry of prosthetic reconstruction, which is currently a key issue.

Brånemark's implantology has yielded exceptional results, but science is a process of both evolution and progress and this process must necessarily begin with and be led by biology. It will be necessary to go beyond the beautiful clinical images and identify clinical and radiological assessment criteria for peri-implant health and its long-term sustainability, which is to define with clear parameters what we might call the "supracrestal integration" of the implant-prosthetic unit. I have enthusiastically joined the Leader Medica research group and its management team for this "Biologically Oriented" vision and I am confident that, with these premises, interesting results will emerge and attract international attention.

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Modern medicine is moving toward minimally invasive and personalized treatments, and dentistry, particularly implantology, is following this trend. The goal is to reduce surgical trauma and optimize tissue healing. The concept of minimally invasiveness has influenced the surgical phase, with computer-guided planning, new surface topographies, and modern macrogeometries that allow for faster healing with less inflammation.

For a long time, prosthesis was considered a secondary step in implantology, but it must be treated as an integral part of the surgical process. The integration of a convergent intramucosal collar improves the response of the peri-implant tissues and contributes to the maintenance of the implant-prosthetic compartment. The primary advantage is the elimination of bone and connective tissue contamination, typical of two-stage implants. The micropumping phenomena occur at the IAJ (implant abutment junction), pushing bacteria into the peri-implant sulcus, increasing the risk of peri-implantitis. The second advantage is the convergence of the neck and the FMC (Fibroblast Micro Channel) surface, which promotes the formation of a connective tissue labrum that prevents apical bacterial proliferation.

The B.O.P.T. technique also allows for the physiological restoration of the anatomy and aesthetics, achieving natural and functional results.*

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Modern implantology increasingly relies on the integration of reconstructive surgical techniques and soft tissue management with dental implants designed to respect each patient's unique tissue healing patterns

Prosthetically driven and computer-assisted protocols for implant placement and the reconstruction of deficient peri-implant tissues (CGBR method) are supported by strong clinical and scientific evidence. Unfortunately, these are often overlooked in favor of outdated approaches that treat surgical and prosthetic phases as separate therapeutic steps. Implant macrostructures and connection designs play a crucial role in achieving and maintaining optimal prosthetic outcomes.

Convergent collar morphologies promote marginal soft tissue healing and enhance long-term tissue stability—especially when combined with B.O.P.T.® prosthetic protocols.

This approach—focused on tissue respect and, where needed, the aesthetically and prosthetically guided reconstruction of missing volumes—represents both the present and the future of modern implantology, grounded in robust clinical and scientific evidence.

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The advancement of digital resources has caused the evolution of implant surgery, which has now been using more precise and reproducible techniques compared to traditional freehand methods. Computer-assisted implant surgery (CAIS) systems enable the detailed visualization of soft and hard tissue, integrating prosthetic restorations and enabling prosthetically guided planning that respects biological, aesthetic, and functional principles.

The Biological Oriented Preparation Technique (B.O.P.T.[®]), developed by Dr. Ignazio Loi, has revolutionized implant treatment by combining tissue modulation with virtual implant planning, gaining increasing relevance in oral implantology. This technique optimizes the final result, increasing predictability and precision in surgical planning. Guided surgery, therefore, becomes a fundamental tool for clinicians, facilitating the integration of restorative factors into virtual surgical planning and improving treatment reliability.

At the University of Valencia, we have spent years perfecting guided surgery protocols, combining our experience in peri-implant tissue management with virtual prosthetic design to improve the long-term stability of tissues. With the support of Leader Medica, we are developing advanced technological solutions to optimize the efficiency and quality of implant treatments with a particular focus on the B.O.P.T.[®] technique and prosthetically guided surgery.

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The Biologically Oriented Preparation Technique (B.O.P.T.®) provides for a new concept in dental prosthetics, not only restoring function and aesthetics, but also shaping and stabilizing soft tissue over the long term.

Our team at the University of Valencia has spent over a decade studying and perfecting this technique, both on teeth and implants, in collaboration with Dr. Ignazio Loi, creator of B.O.P.T.®. This minimally invasive, predictable, and scientifically based approach has led to significant advances in implant treatment.

Intensive research has focused on the evolution of periodontal and peri-implant tissues, marginal stability, and patient satisfaction, and it has confirmed that B.O.P.T.® is a cross-cutting therapeutic technique to be used in restorative dentistry and implantology.

With the support of Leader Medica, we are now embarking on a new, enthusiastic phase, developing simplified and systematic clinical protocols to extend the B.O.P.T.® technique to modern restorative dentistry. We believe that the future of dental and implant prosthetics lies in techniques that respect and enhance biology.

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B.O.P.T.[®] and Loi; Loi and B.O.P.T.[®].

Two words that summarize more than a decade of my professional life. They represent a way of understanding that prosthetics, both tooth- and implant-supported, transcend preparation or restoration: a biological philosophy that redefines the relationship between tissue, prosthesis, morphology, and technique.

From the very beginning, together with Dr. Ignazio Loi, we understood that the goal was not simply to define a finish line, but to free ourselves from it. That freedom has opened a new horizon of biological design, where each restoration can engage in dialogue with the tissue rather than impose itself upon it.

The real learning has been to redirect the clinical gaze toward where stability truly occurs: the pericoronaral mucosal seal. In both teeth and implants, that delicate ring of tissue is the true guardian of the underlying bone.

Today, after more than ten years of applying, understanding, and teaching this therapeutic philosophy, I consider it a true privilege to be part of the research group promoted by Leader Medica. The most extraordinary thing about B.O.P.T.[®] is not what we already know, but what we still have yet to discover.

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B.O.P.T. is a simplified prosthetic procedure that uses vertical preparation in its initial stages.*

The technique was developed in the 1990's through the observation of clinical cases and related events that affected, to varying degrees, the horizontal space and the marginal component of the soft tissue.

The observation, combined with the improvement of the technique, confirmed the initial impressions that the manipulation of these spaces led to a volumetric, and therefore BIOLOGICAL, change in the soft tissue, promoting their medium-to long-term stability.

The two main components of the procedure are the creation of SPACE, through preparation, and SHAPE, through the restoration and its angular component (emergence angle).

Through the use of these two components, we can modulate this space horizontally and sagittally, allowing the soft tissues to grow, occupy the space and increase their thickness. This new HORIZONTAL BIOLOGICAL AMPLITUDE we have named the BIOLOGICAL POLYGON.

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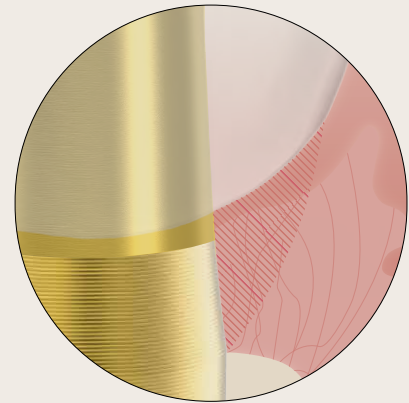
/ THE BIOLOGICAL POLYGON

The BIOLOGICAL POLYGON is determined by the crown's emergence profile, both horizontal and subsequently vertical, that will occupy these spaces, allowing the tissues to remain stable over time. These concepts were then transferred to implantology.

Initially, the construction of converging vertical abutments allowed for increased space, thus also creating a BIOLOGICAL IMPLANT POLYGON for the growth of peri-implant soft tissue and a new implant-prosthetic relationship.

The final evolution of these concepts was the creation of a convergent implant with a neck of different heights. This not only increased tissue volume, but also allowed the implant-abutment junction to be shifted coronally, distancing bacteria from the bone crest.

The IL® implant is the result of this long journey that has always been focused on respecting the tissue BIOLOGY. Numerous published scientific studies indicate that we are on the right path.

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Morphology in comparison

Differences between traditional implant and convergent-profile implant.

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/ MAXIMUM SURGICAL FLEXIBILITY

IL® and IL WIDE® implants are distinguished by the unique shape of their implant body and threads, designed to ensure excellent surgical and prosthetic performance in every clinical situation. Thanks to their geometry, the implants are suitable for placement in all bone types according to the Mish classification: D1, D2, D3, and D4, offering the clinician exceptional operational versatility.

Another distinctive feature of IL® and IL WIDE® implants is that they come in four different intramucosal neck heights: 0.9 mm, 1.8 mm, 2.8 mm and 3.8 mm.

This feature allows the clinician to adapt the implant to both different clinical conditions and different soft tissue thicknesses, choosing the most suitable configuration that best suits the soft tissue morphology. The major advantage is that, regardless of neck height, the same surgical kit and prosthetic components can be used, simplifying their management. This makes IL® and IL WIDE® implants the ideal solution for those seeking efficiency, operational simplicity, and clinical flexibility, without compromising on quality.

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/ TRUST THE ORIGINAL

B.O.P.T.® is an internationally registered trademark since 2014 by Ignazio Loi, exclusively licensed to Leader Medica.

Leader Medica adopts and promotes a protocol aligned with the B.O.P.T.® (Biologically Oriented Preparation Technique) philosophy, which has revolutionized prosthetic management by shifting the focus from a finishing line preparation to a vertical preparation without a margin, achieving thick, stable, and healthy gums. The experience of Dr. Ignazio Loi has allowed these principles to be successfully extended to implantology.

All the devices in this catalog are manufactured by following and promoting the principles of the B.O.P.T.® technique. With the same philosophy, Leader Medica, together with Dr. Ignazio Loi, organizes exclusive courses for professional training, offering original and official contents: the only ones authentically attributable to B.O.P.T.®.

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*A registered trademark since 2014 by Ignazio Loi,
exclusively licensed to Leader Medica*

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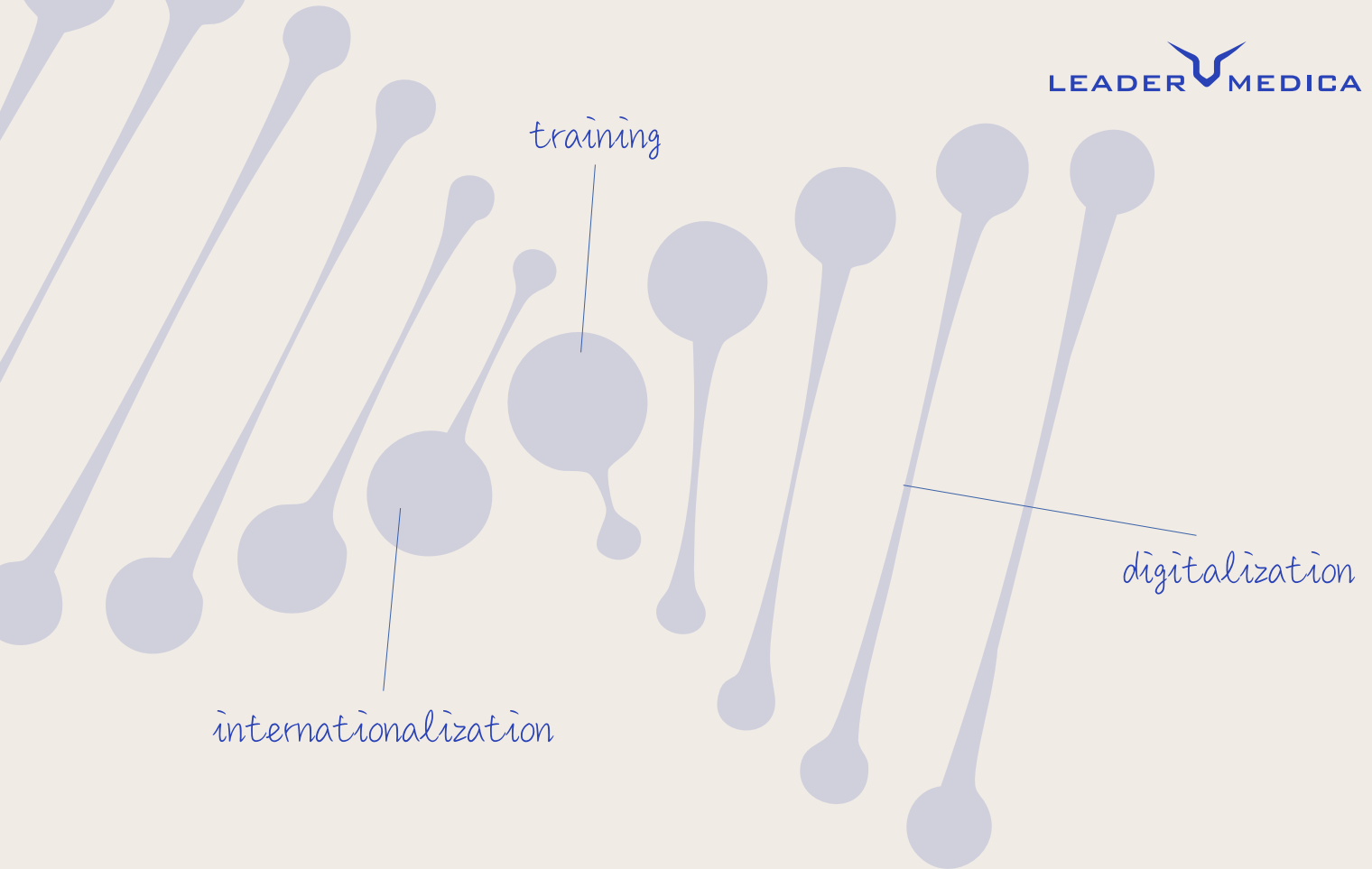
The DNA of Evolution

Leader Medica pursues a Project based on the following core concepts:

- **Maximizing *Quality*** in every phase of the dental implant devices' designing and manufacturing;
- **Training** for anyone who interacts with the company (i.e. internal collaborators, sales associates, dentists, dental technicians, etc.);
- **Digitalization** of both the processes involved in the production of medical devices and the processes relevant to the use of the devices (e.g., guided surgery systems);
- **Generating and/or strengthening a close *Relationship of Trust*** between the Company and the Client with the aim of creating the foundation for a successful and synergistic business and professional development, characterized by valuable training and digital content.

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Device quality, training, digitalization, and internationalization are the key concepts of the evolution that sets us apart and defines us for the future.

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Leader Medica has been designing and manufacturing high-precision implant systems and prosthetic components since 1994, overseeing every phase in-house, from the research to the production as well as the clinical and training support. Innovation is part of our DNA: we study, discover, and create solutions that simplify and improve therapeutic efficacy. We are a dynamic, constantly evolving company focused on international trade and boasting a network of outstanding scientific collaborations.

LEADER MEDICA consolidates its growth

In 2024, Leader Medica acquired 3P Implafavourite, a company that has distinguished itself in the market with the 3P Smart Devices brand and as the creator of the CONO-IN® implant line featuring a conical Morse-taper connection, a true expression of technological and design excellence. Through this acquisition, a second manufacturing facility is established in Pinerolo (Turin, Italy), strengthening the Group's industrial capacity and strategic footprint.

LEADER MEDICA worldwide

2

2 SEDI PRODUTTIVE IN ITALIA
2 PRODUCTION SITES IN ITALY

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115+

COLLABORATORI COMMERCIALI
BUSINESS PARTNERS

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Leader Medica develops excellent training programs, tailor-made for dentists seeking concrete and specialized growth. We offer custom training with a focus on convergent implantology, everything related to the digitalization of dental implantology procedures and protocols, including guided surgery and conometric systems. The training programs are designed with the involvement of the attendants, ranging from introductory courses for beginners to more advanced courses, also including practical applications.

Leader Medica stands out for its decisive international focus, consolidated by distribution partnerships in several countries around the world and collaborations with prestigious academic institutions such as the University of Padua, Tuscan Stomatological Institute, Magna Graecia University, University of Genoa, University of Seville, University of Valencia.

Training and Digitalization, with a multifaceted international vision and culture, are an **integral part of our Project** that is the **DNA of the evolution** that we wish to carry out together with you.

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*From the origin of convergence,
IL[®] original convergent*

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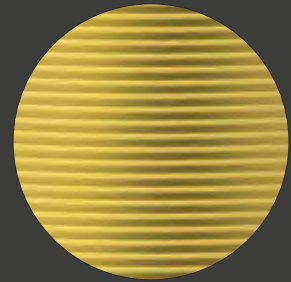
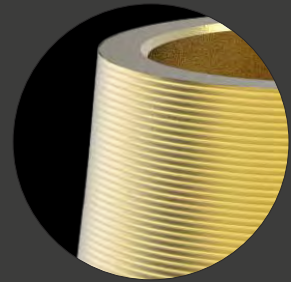
FMC NECK:
FIBROBLAST MICRO CHANNEL

STRAIGHT AHEAD:
NECK SHAPE

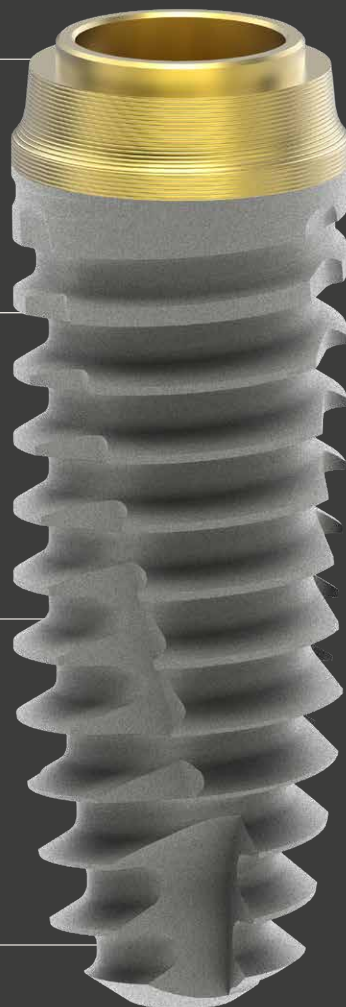
THREADS: ATC ASYMMETRICAL
TRAPEZOIDAL COIL

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*From the origin of convergence,
the original convergent IL[®]*

01

COLLO CONVERGENTE HPA | HPA NECK HIGHT PERFORMANCE ANODIZING

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The yellow-gold anodic coloring of the neck improves the biomimetic capabilities at the intramucosal level.

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The implant body is slightly tapered to support and improve primary stability.

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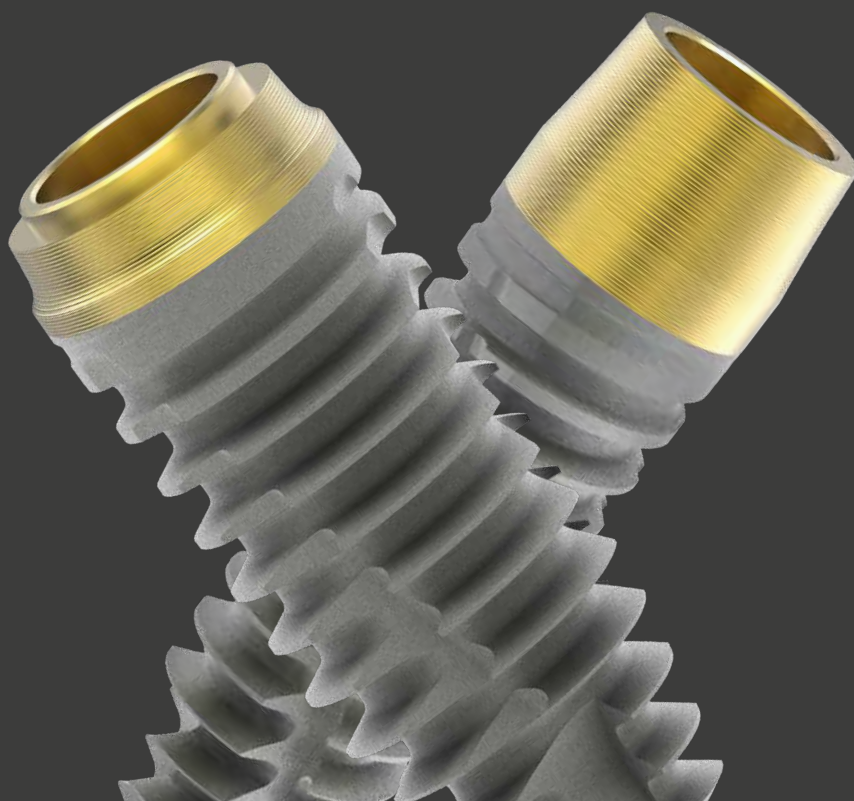
The external dimension of the ridges is tapered with a specifically designed taper angle to facilitate insertion during osteotomy.

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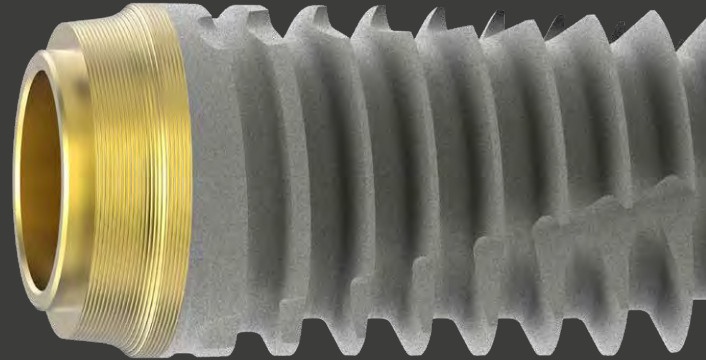
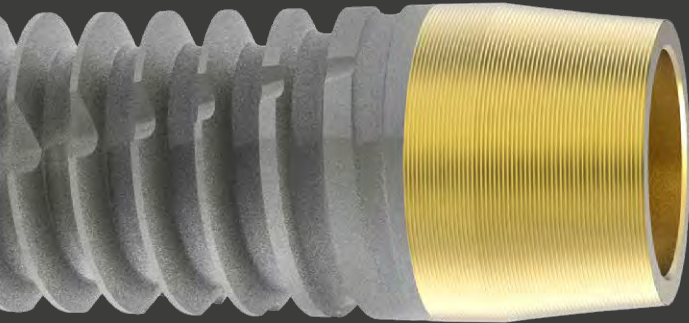
They facilitate the initial indentation of the thread, ensuring immediate correct engagement of the implant in the bone and therefore its sinking into it.

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*The convergent neck implant is available with two internal hex prosthetic connections:
 IL with Ø 3.40 mm, a standard connection,
 IL WIDE® with Ø 3.80 mm, specifically designed for posterior restorations.*

The Ø 3.80 mm connection features an internal hex and a 0.32 mm external collar, which increases the mechanical stability of the connection and reduces the stress on the prosthetic screw.

IL WIDE® also favors a more harmonious emergence profile because the connection platform is wider in the cases of posterior teeth prosthesis and that improves aesthetics and functionality.

The use of only two connections allows for high prosthetic versatility, simplifying clinical management and ensuring reliable solutions.

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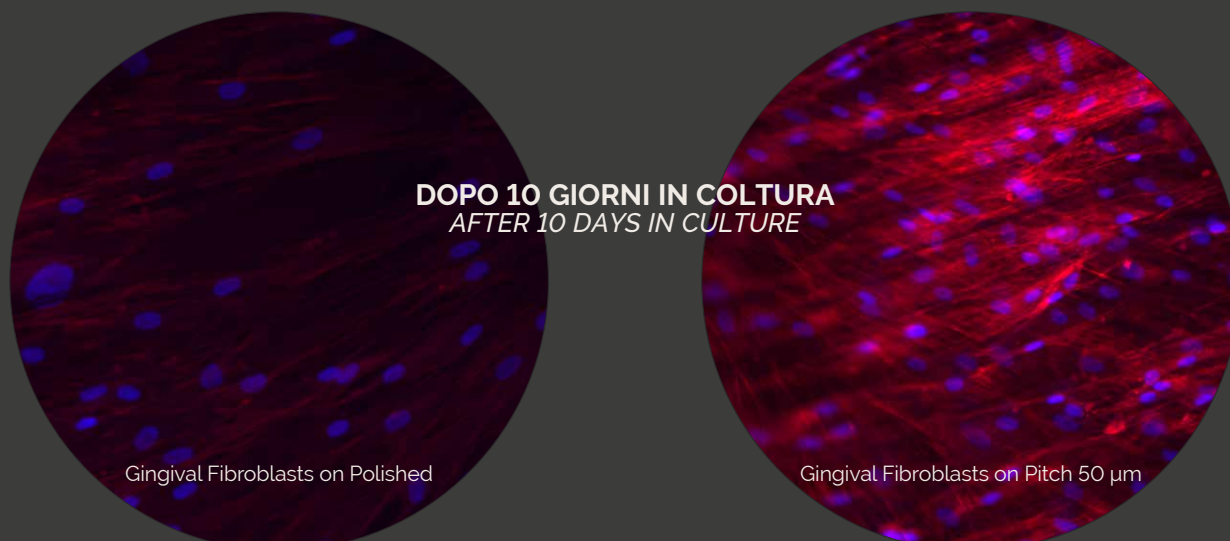
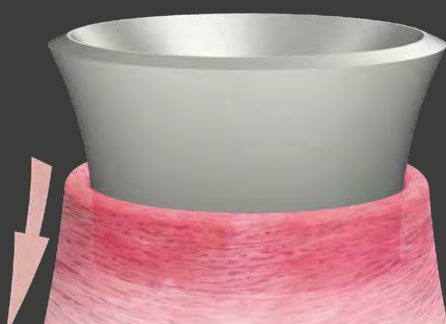
FMC® SURFACE - FOR A LONG-LASTING BIOLOGICAL SEAL

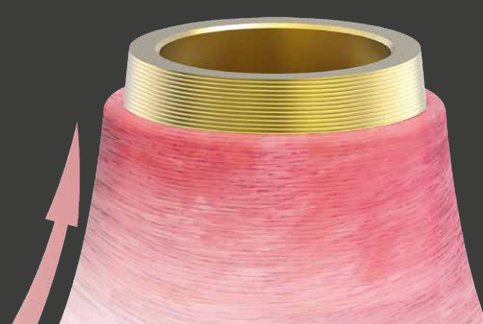
Thanks to advanced fluorescence microscopy techniques and biochemical analysis, it was possible to observe how gingival fibroblasts, endothelial cells, and mesenchymal cells of dental origin interact selectively with different implant surfaces. The results demonstrate that micro-grooved surfaces 50 μm , such as the FMC® (Fibroblast Micro Channel) surface of the neck of IL® and IL WIDE® implants, promote better cell adhesion and activation, increasing collagen deposition, a key element for soft tissue stability. This contributes to the formation of a stable biological seal, resulting in faster healing, and it ensures safer, longer-lasting, and more aesthetically natural results for patients. Optical microscopy analysis demonstrated a highly positive effect on the convergent micro-grooved surface with greater density and improved interweaving of the collagen bundles. Pending research has also highlighted the ability to stimulate the transition from fibroblasts to myofibroblasts, likely enhancing the stability of the connective tissue seal.

The study was conducted under the aegis of Professor Di Vito, in full compliance with the most rigorous methodological standards.

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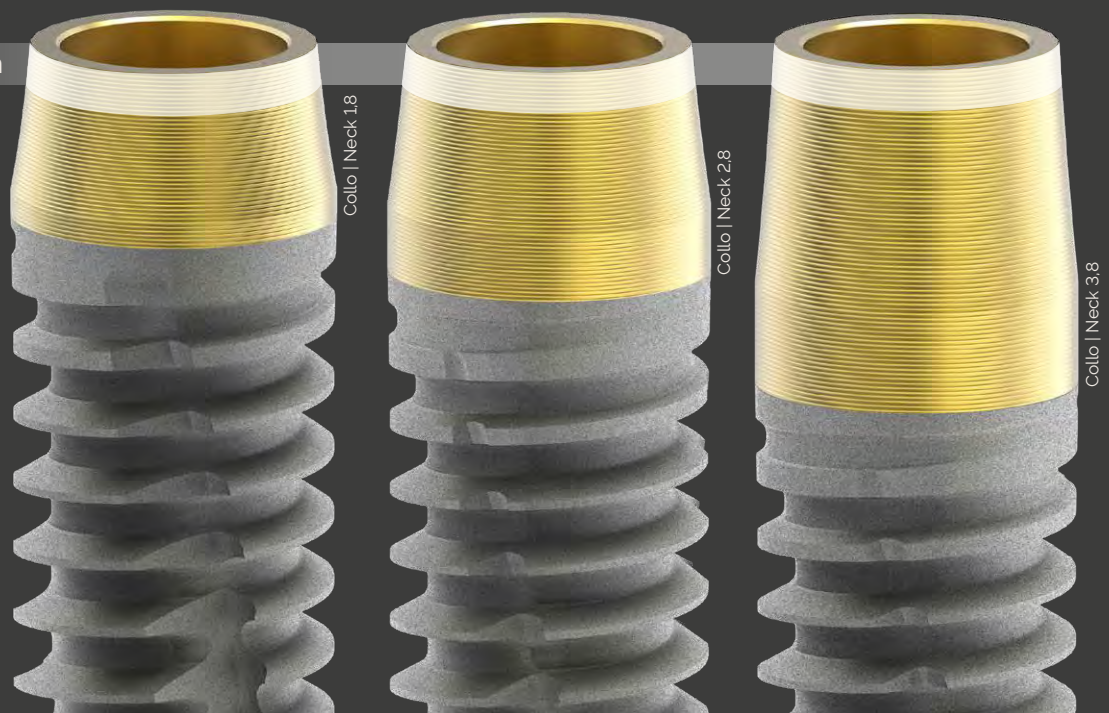
The specific conical shape of the first 0.5 mm of the implant neck ensures perfect screwing of the Biologic abutment onto the intramucosal neck of the implant, expanding the contact surface. This, especially in single elements, makes the connection much more resistant and stable. In cases with thinner soft tissue, it allows for the insertion of implants using the B.O.P.T.® technique to promote gum's horizontal thickening. This allows the use of the same prosthetic components across the entire IL® implant line.

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**REDUCED RISK OF
CONTAMINATION WITH
FOREIGN MICROPARTICLES**

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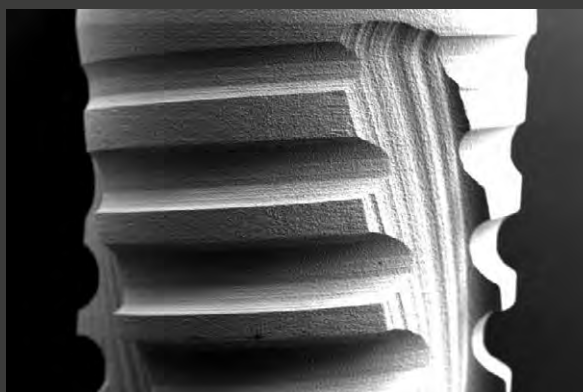
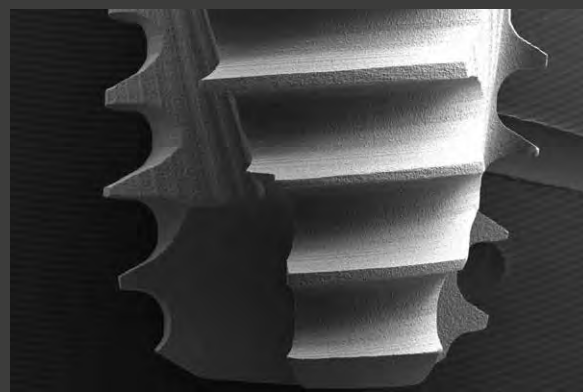
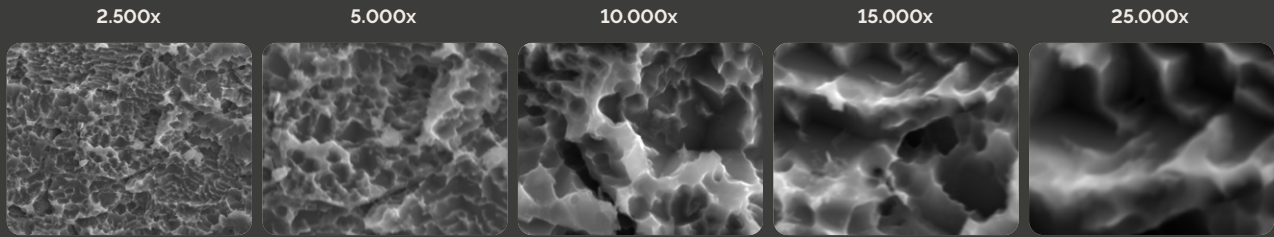
The DAS surface (present on the body of all Leader Medica implants) is a double-acidified surface. This surface is obtained using a chemical method that involves the use of acid mixtures to roughen the titanium surface by subtraction. The first phase involves a surface activation treatment (i.e. the creation of a slight macro-roughness) using hydrofluoric acid. After a thorough washing with demineralized water, the actual etching responsible for the micro-roughness occurs using a mixture of sulfuric and hydrochloric acids. After a further wash, the passivation process takes place: it's a deep cleaning of the surface to restore the passive layer of the material through the controlled supply of oxygen that promotes the formation of titanium oxide. This type of surface allows for an arithmetic mean roughness (Ra) of 1.35 μm . These roughness values are comparable to those expressed by SLA-type surfaces (Sandblasted Largegrid and Acid-etched) and are in perfect agreement with the literature. As demonstrated by Wennerberg A., the long-term results between these two moderately rough surfaces are superimposable. The advantage of the DAS surface (which is much more costly than the SLA one) is that it prevents the deposition of blasting particles on the implant surface, which often results in potentially harmful aluminum oxide residues.

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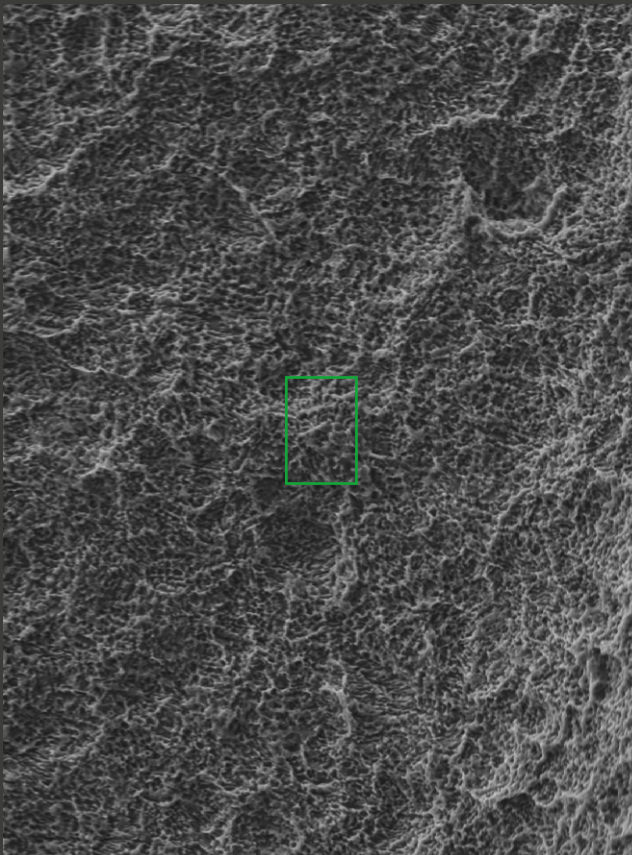
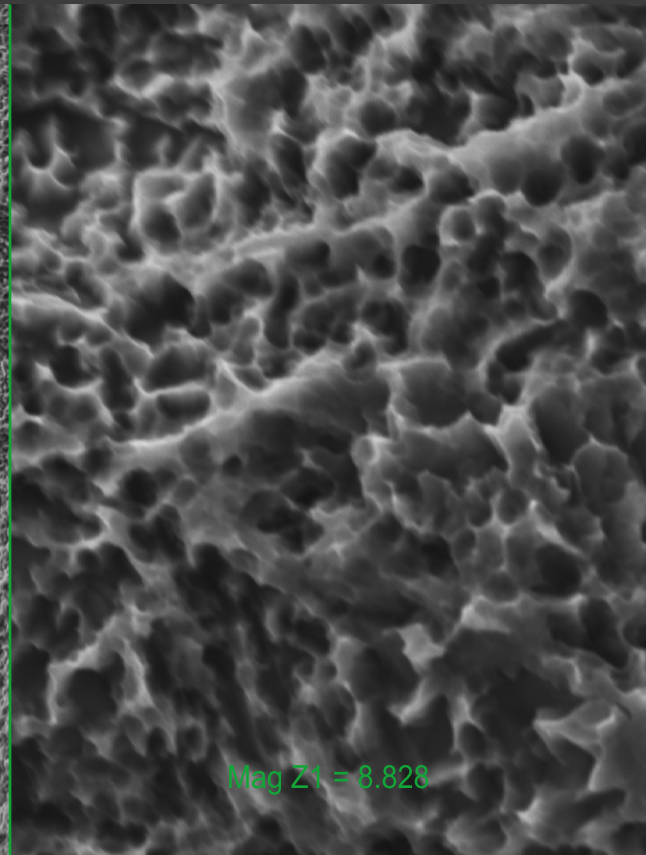
We internally certify surface roughness and cleanliness by using scanning electron microscopy (SEM) and X-rays.

The protocol applied to each batch of implants was defined with the contribution of Professor Luigi Canullo, recognized as one of the leading experts in the field of implant surfaces and related manufacturing and cleaning processes.

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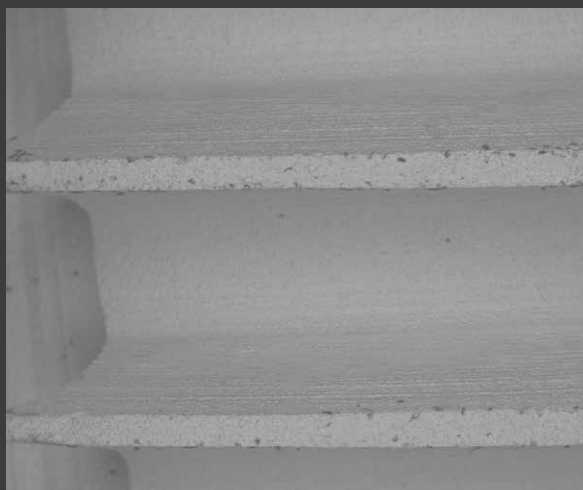
Industrial processes (turning, surface treatment, packaging) always produce processing residues that deposit on the surface. Although decontamination phases involving washing and ultrasonic baths are included in the various steps of the manufacturing process, microscopic examinations often reveal carbon residues or acid residues. The effect of these contaminants undoubtedly decreases osteoblastic cell adhesion and simultaneously produces an inflammatory response.

The most effective strategy for removing them is to use a biophysical process called plasma. This technology uses an ionized gas or an ionized gas mixture to remove any contaminants from the implant surface.

The plasma used for Leader Medica implants, through specific settings and argon-oxygen gas mixtures, effectively removes organic and inorganic contaminants, making the surface ready for interaction with bone cells.

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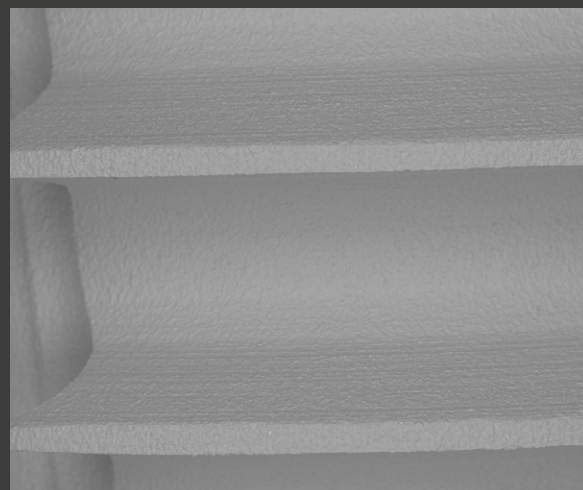
IMPLANT BEFORE TREATMENT

SEM | Mag- 57 X Signal A- BSD1

IMPIANTO DECONTAMINATO

DECONTAMINATED IMPLANT

SEM | Mag- 57 X Signal A- BSD1

SPIRE ATC - Asymmetrical Trapezoidal Coil

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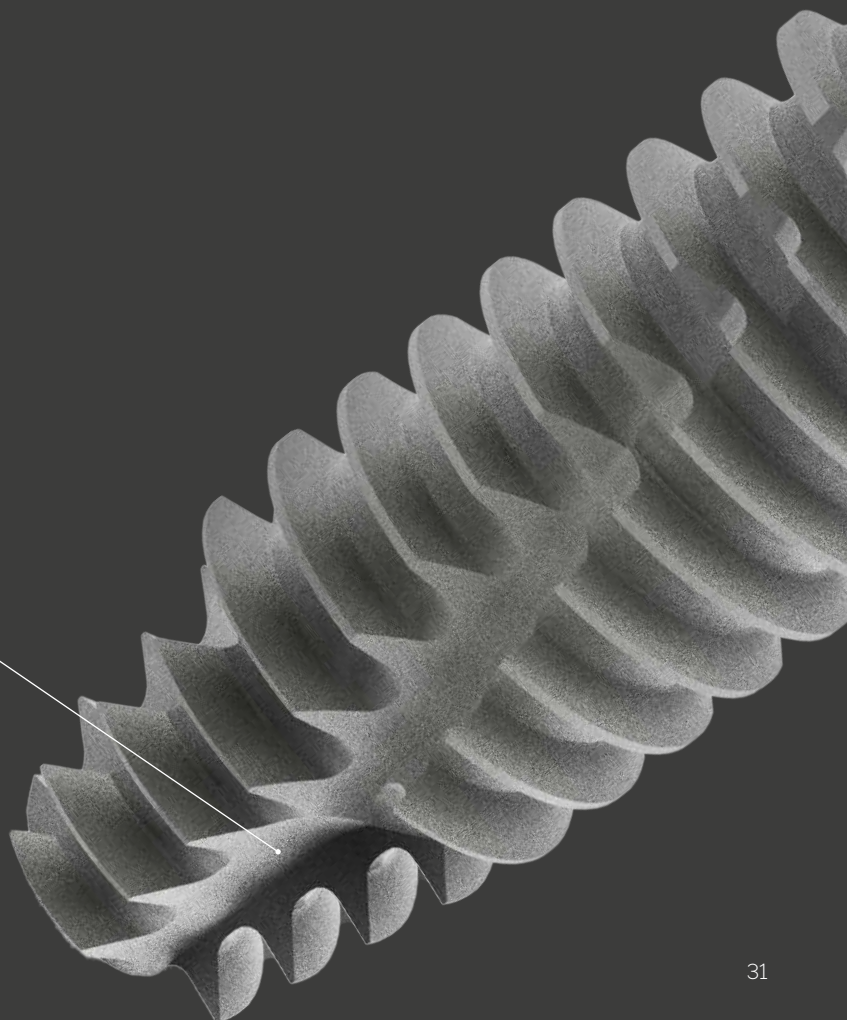
The asymmetrical trapezoidal coils improve primary stability. This shape ensures the best distribution of masticatory loads in the bone. Their increasing depth in the coronal-apical direction allows for excellent implant stability, even in predominantly medullary bone.

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The blades (Bodyblade®) are used to tap the socket, reducing the insertion torque. In post-extraction implantology, the blades work like hooks, clinging to the (asymmetric) alveolar wall, allowing the surgeon to maintain the set direction without deviations due to the reduced resistance of the bone where it is less compact and stable or towards the alveolar cavity.

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Certification of the Protocol for Surface Treatments and Processes

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WHAT IS THE LEADER MEDICA CERTIFIED SURFACE?

The LEADER MEDICA CERTIFIED PROTOCOL concerns the processes and treatments performed on the surfaces of Leader Medica dental implants, as shown in the diagram on the following page.

For the purpose of implant osseointegration, the considered implant surface parameters are, on the one hand, ROUGHNESS, and on the other, CLEANLINESS and PURITY.

Roughness refers to the deviations of the peaks and valleys of the surface asperities calculated along a linear path traced on the surface;

Cleanliness and Purity refer to the absence of localized accumulations of contaminants on the implant surface.

CREATED BY

The LEADER MEDICA CERTIFIED PROTOCOL was developed, refined, tested and validated by the engineers, chemists and technicians of Leader Medica, together with Professor Luigi Canullo (Unicamillus University, Rome).

WHY?

The LEADER MEDICA CERTIFIED PROTOCOL was developed to enhance the quality of surface Roughness, Purity, and Cleanliness, ensuring the predictability and repeatability of these parameters.

IS IT TRUE?

Leader Medica have a latest-generation SEM (Scanning Electron Microscope) equipped with EDS (Energy Dispersive Spectroscopy), used to verify and capture images of the surface Roughness and Cleanliness of dental implants treated according to the LEADER MEDICA CERTIFIED PROTOCOL.

The Roughness images (with 8000x magnification) and Surface Cleanliness images (via energy-dispersive spectrometry) are produced by Leader Medica and always provided to the clinician together with the dental implant belonging to the batch verified by SEM-EDS.

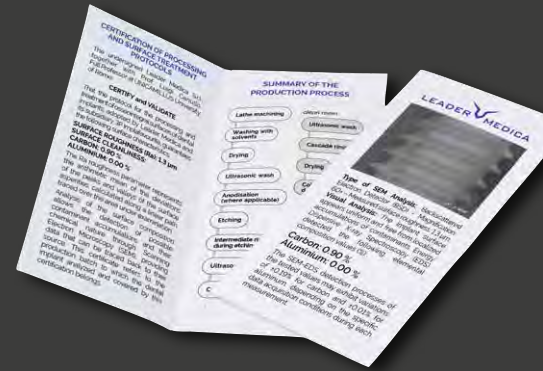
FINAL RESULTS

The results obtained through the processes and treatments of the LEADER MEDICA CERTIFIED PROTOCOL, in terms of Roughness, Purity, and Surface Cleanliness, are considered EXCELLENT:

1. Surface roughness (Ra) measured at 1.3 µm;
2. A characterization of absolute purity and surface cleanliness: (Aluminum: 0.00% and Carbon: 0.90%, based on a verified sample of 5,000 dental implants)

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02

LAVAGGIO CON SOLVENTI
Washing with solvents

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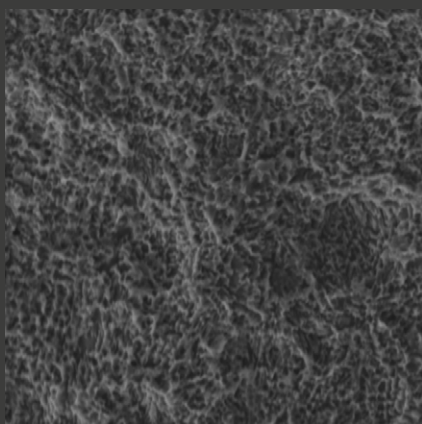
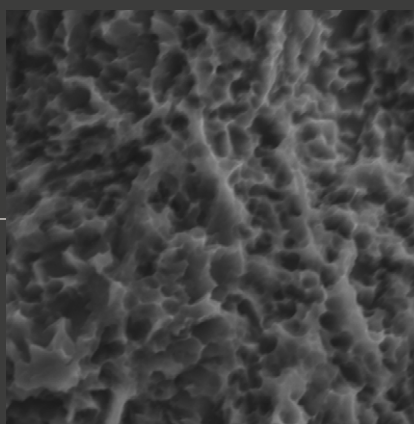
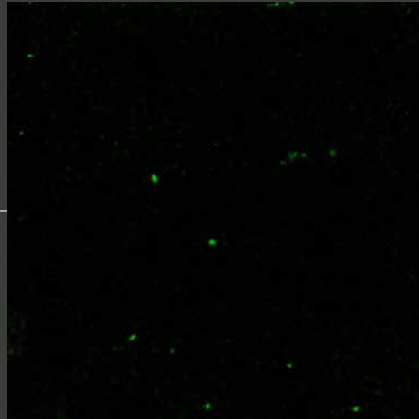
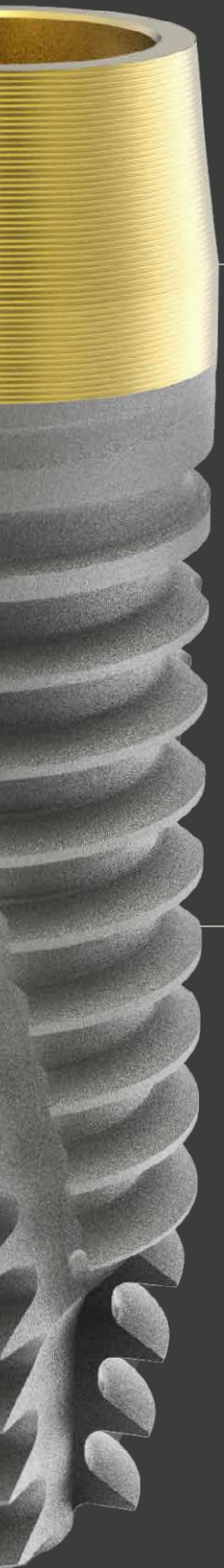
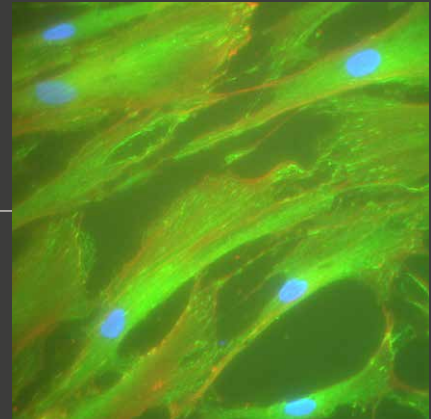
PULIZIA SUPERFICIE
SURFACE CLEANING

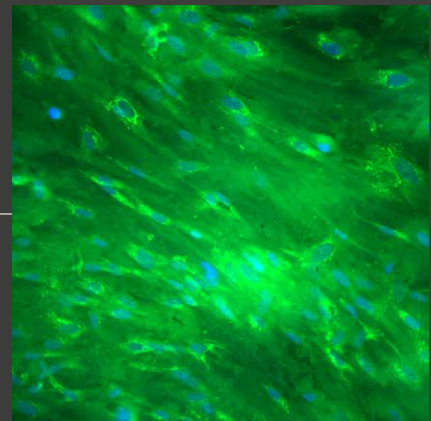
SEM | EDS Mag- 60 C 0,08%

Gingival Fibroblasts on Pitch 50 μ m. Vinculin (green), F-actin (red), nuclei (blue), 40X

ADESIONE CELLULARE
CELLULAR ADHESION

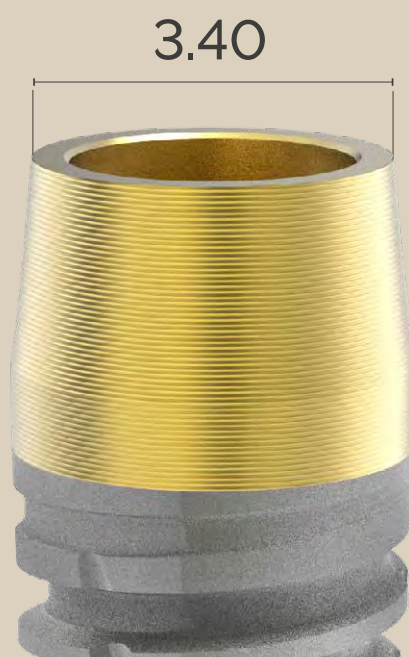
Dental mesenchymal cells. Collagen I (green) and nuclei (blue), 20X

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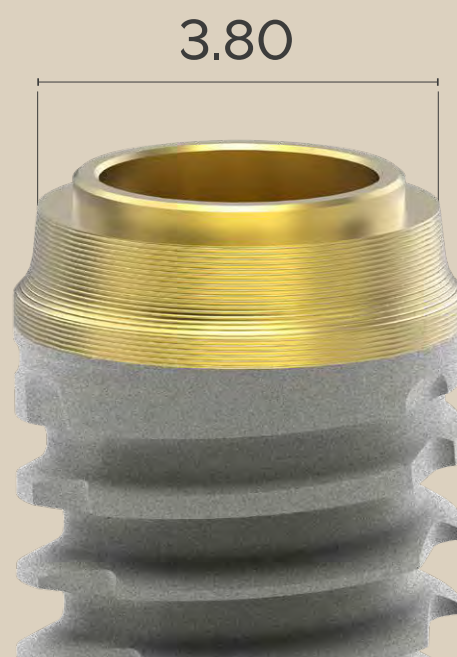
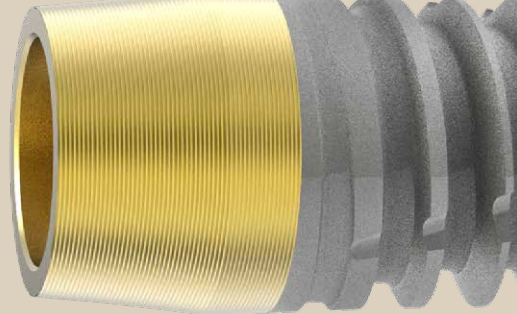
IMPLANT LINE IL®, PROSTHETIC PLATFORM CONNECTION 3.40 mm

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IMPLANT LINE IL WIDE®, PROSTHETIC PLATFORM CONNECTION 3.80 mm

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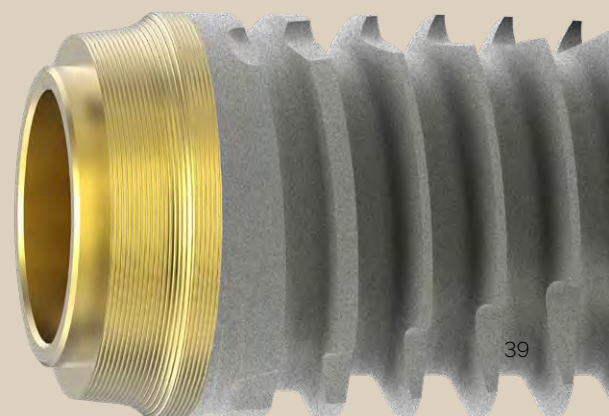
Surgical protocol, insertion sequence

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The principles of B.O.P.T. have been successfully extended to implantology, creating a specific surgical sequence for the IL* and IL WIDE* implant lines. IL* implant placement guidelines take into account soft tissue thickness and the anatomical conformation of the bone crest. This integrated approach allows for a prosthetically guided relationship in implant placement. The formation of circular collagen fibers in the microgrooves of the convergent neck, along with the formation of a short epithelial attachment around the prosthetic device, guarantee the long-term stability of the IL* implant.*

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For illustrative purposes, the surgical sequence for an implant with the following characteristics is shown below:

**Ø 3.80 MM, H 10.8 mm, NECK 2.8 mm - STOP USED
FOR THE 10 mm HIGH IMPLANT.**

	PRIMA FRESA FIRST DRILL	SECONDA FRESA SECOND DRILL	TERZA FRESA THIRD DRILL	QUARTA FRESA FOUR DRILL	QUINTA FRESA FIFTH DRILL
Ø EFFETTIVO FRESA ACTUAL DRILL	LANCEOLATA LANCEOLATE	2.0	2.5	2.8	3.0
Ø IMPIANTO IMPLANT	-	-	-	3.4	3.8

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IMPLANT INSERTION

The IL[®] implant is a tissue-level implant and therefore must be placed with the convergent neck portion outside the bone. This maximizes the advantages of having the abutment/prosthesis connection far from the bone and the collagen fibers can create a biological seal around the implant neck.

The clinician decides on the placement based on the quantity and quality of the soft tissue, taking care to leave the connection slightly below the gingival margin.

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NECK CHOICE

The bone crest and gingival thickness can vary when choosing between intramucosal or tissue-level insertion that is why we have introduced four neck heights: 0.9, 1.8, 2.8, and 3.8, to adapt to any situation.

The implant can be placed in a tissue-level position with minimal controlled soft tissue invasion, or in a slightly deeper intramucosal position. The decision will be made by the clinician during the case planning phase, also taking into account the prosthesis with which the patient will be rehabilitated (i.e. cemented, screw-retained, single or multiple implants). Correct implant placement should take into account the CEJ of the adjacent tooth, if it is present, aligning the connection with this line. In the case of more extensive edentulism, the location of the new crown must be considered so as to respect the concept of biological width (the minimum mandatory requirement is approximately 2 mm).

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USE OF THE STOP AND THE COMPENSATION SPACE

The use of the stop is designed to allow the first 0.8 mm of the neck to sink into the bone.

This surface, called the compensation space, was designed to adapt to and compensate for the inevitable post-surgical bone remodeling.

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THE IL WIDE[®] IMPLANT

For greater stability and mechanical strength in the posterior quadrants, we recommend using the IL WIDE[®] implant with the same insertion procedures as the IL[®] implant.

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THE COMPENSATION SPACE

When measuring bone availability on radiographic images, consider that the drill action includes 0.8 mm, which must be added to the size of the chosen implant. The 0.8 mm refers to the compensation space provided by the protocol.

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TABLE SUMMARISING DRILL DIAMETER AND LENGTH

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To accurately determine the drilling depth, it is recommended to consider the length of the drill bit, which varies for each diameter.

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In keeping with the close synergy between the B.O.P.T.[®] philosophy and the convergent implantology of the IL[®] line, Leader Medica has developed a series of dedicated B.O.P.T.[®] preparation kits. These kits were designed in direct collaboration with Professor Ignazio Loi and manufactured according to his specific clinical indications and operating protocols with the aim of offering perfectly integrated and biologically oriented surgical and prosthetic solutions. Each drill has a specific use in the various stages of preparation up to the creation of the prosthetic abutment.

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The vertical preparation is created, using dedicated diamond drills of various shapes and grit, depending on the step of the preparation they are going to be used for.

PREPARATION AND FINISHING DRILLS FOR TEMPORARY BRIDGES

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The sequential use of drills, discs, rubber, and brush allows for specific use in the various steps until the final, finished and polished temporary restoration is created. Each drill has a specific function in each step.

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14

STOP PER FRESA Ø 2,0 E Ø 2,5
STOP FOR DRILL Ø 2.0 AND Ø 2.5

- ILSTOP20-2504*
- ILSTOP20-2506
- ILSTOP20-2585
- ILSTOP20-2510
- ILSTOP20-2511
- ILSTOP20-2513
- ILSTOP20-2515

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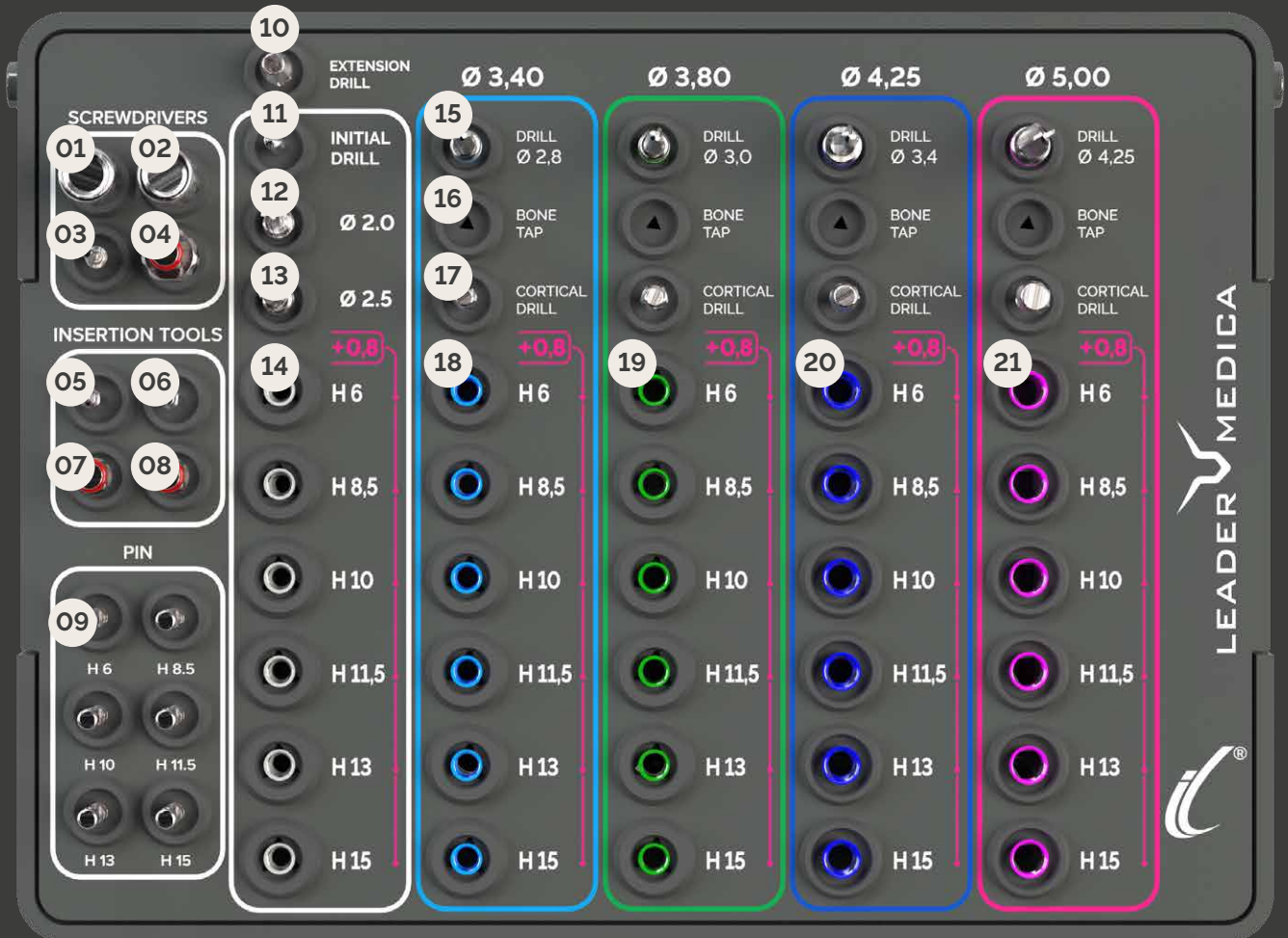
From analog to digital...

Precision, clinical predictability and aesthetics are at the core of our digital workflow. Each step of the process is carefully designed to effectively meet the needs of patients who require rapid and minimally invasive solutions with short recovery times, as well as professionals who seek maximum operational safety and reliable results, both functional and aesthetic.

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“ *The meeting point between surgeon,
dental technician, and patient* ”

Thanks to virtual planning integrated with 3D intraoral scanning, it is possible to acquire extremely accurate data, faithful to the patient's anatomy, and develop a personalized prosthetic design. Prosthetically guided surgery enables more controlled, predictable, and time-optimized procedures, improving the clinical experience and patient comfort.

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Data acquisition via intraoral scanner STL format and CT or CBCT in DICOM format.

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Implant and prosthetic design using 3PGuide software.

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3D design and production of the surgical guide, anatomical model, and temporary prosthesis.

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Prosthetically Guided Implant Surgery.

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With the 3P Guide software from Leader Medica, we provide you with a powerful tool for a comprehensive diagnosis and extremely accurate surgical and prosthetic planning, allowing you to complete the entire workflow.

Flexible: Compatible with all DICOM formats, it easily adapts to any radiological examination. It is also compatible with Mac, Windows and iOS.

Integrable: it interfaces with major CAD/CAM systems for fluid and consistent modeling.

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Interconnected modules for a seamless workflow

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Integrated cloud-based file sharing platform

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Intuitive procedures, thanks to automated functions

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Open architecture for maximum flexibility

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Integrated artificial intelligence

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Why choose surgically guided implantology?

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01

SAFETY DURING SURGICAL PROCESSING

The combination of intraoral scanning and 3D radiology, integrated into digital planning, allows for an accurate analysis of sensitive anatomical structures (such as nerves, vessels, and bone structure), ensuring safer, more controlled, and highly predictable procedures.

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02

PREDICTABLE FINAL RESULTS

Digital tools allow for precise visualization of the final prosthetic results already in the pre-operative phase, improving the planning and effectiveness of the entire treatment.

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03

MAXIMUM OPERATIONAL PRECISION

The customized surgical guide, created with a 3D printer, follows the patient's anatomy, ensuring the correct implant placement and reducing the margin for error and post-surgical complications.

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04

SHARED CLOUD SOFTWARE

The integration of digital platforms with cloud access facilitates seamless and immediate communication between the clinician and the laboratory, allowing for real-time project sharing.

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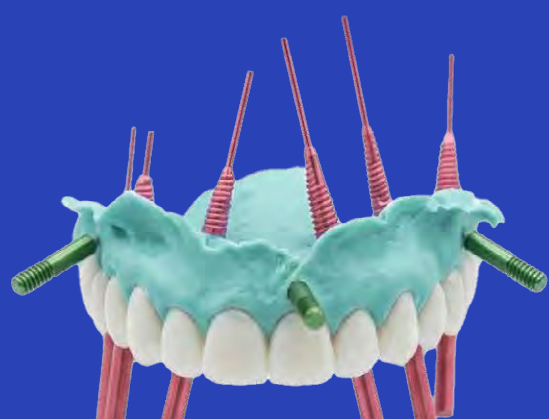
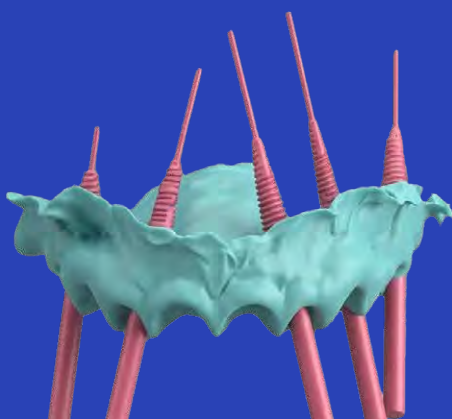
05

DEDICATED TECHNICAL SUPPORT

A team of specialists provides ongoing assistance for the use of digital systems.

Protocollo Chirurgia Protesicamente Guidata

Prosthetically Guided Surgery Protocol

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MEDICAL HISTORY AND DENTAL EXAMINATION

The medical history is used to collect all the patient's data to allow for a correct diagnosis. The dental examination aims to assess the general health of the mouth, preventative care, and assess the operating spaces, as a limited oral opening can impede the correct execution of the technique or limit the passage of drills. An opening of at least 40 mm is recommended.

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DATA ACQUISITION - PARTIAL EDENTULOUS PROTOCOL

FILE ACQUISITION FROM AN INTRAORAL SCANNER.

The digital process involves the use of an intraoral scanner to improve precision and patient comfort. Digital impressions can only be taken in cases of partial edentulism, with the goal of obtaining a suitable oral record of the soft tissues and tooth surfaces.

Scan files must be saved in STL/PLY/OBJ format.

To achieve optimal scan quality, follow the specific protocol for the intraoral scanner used. Scans must be free of holes and well-extended over the fornix and palate in the case of the upper arch.

FILE ACQUISITION FROM A 3D SCANNER

Starting with an analog workflow, the initial plaster model is taken using high-precision material, ensuring full extension, including the palate and retro-molar tuberosities.

The resulting plaster model is then scanned using a benchtop scanner, with careful attention to proper alignment with the indexes and, if applicable, with the prosthesis, reference markers, or diagnostic wax-ups. The entire dataset must be exported properly aligned in .STL format.

DIAGNOSTIC WAX-UP (IF APPLICABLE)

The diagnostic wax-up provides information about the crowns to be restored.

Having a diagnostic wax-up available for the case will simplify planning the aesthetic aspects of the implants and abutments, allowing for correct guided prosthetic planning. It can be created directly in digital format or derived from a scan of an analog wax-up.

The required digital format for the wax-up is STL.

CT SCAN ACQUISITION

The CT scan is performed by acquiring images of the arch, taking care to place special instruments between the two arches to avoid occlusion.

CT OR CBCT - DICOM FILE

The volumetric scan provides information about the bone, teeth, and other anatomical structures.

When scanning a single arch, the maxilla or mandible and part of the opposing bone must be completely visible.

This will make aesthetic assessment easier during planning review. Making the ostiomeatal complex visible will also be useful for evaluating sinus pathologies, if any.

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VOLUMETRIC SCAN SETTINGS:

Scan images must be saved in DICOM format. To achieve optimal image quality, check the device's specific DICOM export protocol. If an export protocol for volumetric scanning is not available, it is recommended to save DICOM images with a recommended axial slice thickness of 0.2 or 0.4 mm.

DATA ACQUISITION - TOTAL EDENTULOUS PROTOCOL WITH DOUBLE SCAN

PREPARATION OF THE SCAN PROSTHESIS

The scan prosthesis is a pre-existing removable prosthesis or a radiolucent duplicate of the diagnostic dental setup with eight markers inserted or glued to the surface (e.g., quartz spheres, Fig. 1-2). The scan prosthesis is used to visualize the desired final prosthetic goal and to design the surgical guide.

1. Verify the pre-existing scan denture. A pre-existing denture can be reused if three conditions are met:
 - a) The dental setup is the final goal;
 - b) There are no metal parts, or the metal parts are removable for scanning;
 - c) The denture fits perfectly or can be modified to perfectly fit in the patient's mouth.

Fig. 2

If the above criteria are not met, a radiolucent duplicate of the diagnostic dental setup will be needed. Specifically, a relined denture is required.

2. Create eight small, evenly distributed surface cavities on the denture (on the vertical, buccolingual, and mesiodistal portions).
3. Use a drop of radiolucent resin and dip a marker into it. Place the resin with the marker in one of the cavities. If necessary, light-cure to reduce the curing time. Repeat this step for all markers.

IMPORTANT: The scanning denture must be radiolucent and free of metal parts.

PREPARING THE OCCLUSAL INDEX

The occlusal index stabilizes and separates the maxilla and mandible during the patient scan. Scans showing maxillary and/or mandibular movement cannot be used to design the guide.

How to prepare an occlusal index:

1. Apply radiolucent material to the patient's occlusal surface.
2. The patient should close their mouth, but not completely. The occlusal index should be 3-5 mm thick; subsequently remove any excess material with a sharp instrument.

Caution! Do not discard the occlusal index after scanning; the occlusal index can be reused to accurately position the guide.

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CT OR CBCT - DICOM FILE

The volumetric scan provides information about the bone, teeth, and other anatomical structures. In cases of complete edentulism, two scans are required, as per the double scan protocol.

How to Perform a Double Scan:

Scan 1

Scan the patient with the occlusal index and the scanning prosthesis with the markers.

If you wish to use an external marker, the marker must be assembled with the scanning prosthesis, always securing it with radiolucent material.

The maxilla or mandible, including the occlusal index and part of the opposing bone, must be completely visible. This will facilitate the aesthetic assessment during planning. Making the ostiomeatal complex visible will also be useful for evaluating sinus pathologies, if any.

Scan 2

Scan the scanning prosthesis alone in the same right/left, upper/lower position used in the first scan.

The support material used for the scanning prosthesis must be more radiolucent than the prosthesis itself. Polyethylene and polyurethane foam are types of materials that can be used when performing a cone beam scan.

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MATCHING

After importing the DICOM of the CT scan into the planning software, and before performing the implant planning, it is necessary to register the STLs according to the type of protocol.

PARTIAL EDENTULOUS PROTOCOL

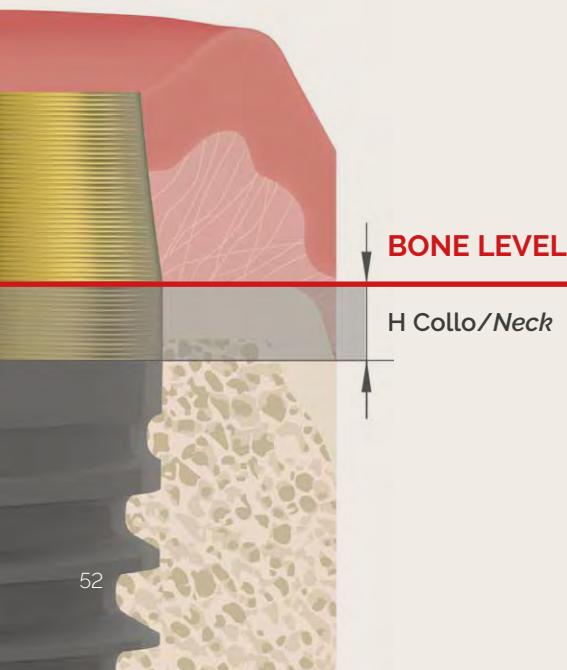
CT - Arch STL (in the case of intraoral scanning or benchtop scanning)

DOUBLE SCAN PROTOCOL

CT (patient + prosthesis + markers) - CT (prosthesis + markers)

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The neck section (H) must be inserted into the bone for implants with a 1.8, 2.8 or 3.8 mm neck (N). Implants with a 0.9 mm neck (N) are bone-level implants.

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The length of the drills and the offset of the cannulae allow for guided implant insertion as outlined in this protocol.

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The design of the drills allows the implant to be inserted into the bone an additional 2 mm compared to the protocol for specific clinical needs.

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The use of the IL® implant with a 3.4 mm diameter must be carefully evaluated by the surgeon. It is indicated exclusively for placement in aesthetic areas, characterized by limited occlusal loading and angulations close to the implant axis.

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3D DESIGN AND PRODUCTION OF THE GUIDE

Based on the planning, the surgical guide is designed. The software allows you to define parameters such as offset, thickness, and undercut management. When designing the guide, it is important to consider the presence and size of the bushings, their spacing, and any support areas that may affect stability or create pre-contact with any existing teeth. In cases of total edentulism, the use of lateral pins and their dedicated drills is mandatory. The surgical guide, designed in the software, is exported in STL format and printed with a 3D printer.

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STANDARD DRILL SEQUENCES AND IMPLANT SURGERY

ANESTHESIA

It is important not to create anesthetic boluses that could prevent the proper seating of the guide.

MUCOTOMY

Performed using a special Tissue Punch after seating the guide.

OPERCULECTOMY

The surgical guide must be removed to facilitate the removal of the mucosal operculum.

CRESTAL PREPARATION AND INITIAL DRILL

It creates the crestal plate for the implant head by removing all interference related to the alveolar ridge. The IL® guided kit includes both a dedicated crestal mill for each type of neck and the initial drill.

DRILLING SEQUENCE

Starting with the first 2x6 mm drill, the subsequent longer drills are guided by both the surgical guide bushings and the 6 mm pilot drill.

The drilling sequence may vary, depending on implant types and the bone quality at the implant site.

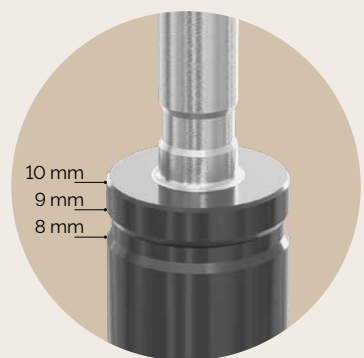
On the left, the different recommended sequences:

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The Crestal Mill has no stop.

It features dual-purpose notch at 8mm and 9 mm. It is possible to reach 10 mm based on the coronal surface of the mill.

IMPLANT PLACEMENT

Special guided mounts screwed into the implant are used to position the implant while maintaining the insertion axis as given by the software design. These mounts also allow you to track the position of the connection system. Different mounts are available, depending on the neck length of the IL® or IL WIDE® implant.

REMOVAL OF MOUNTS AND FIXING SYSTEMS: PROSTHETIC PHASE

Removal of mounts and, if applicable, insert or remove healing screws or first-stage screws.

IMMEDIATE LOADING

The decision to perform immediate loading is at the discretion of the professional performing the procedure.

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Piano inclinabile
Tilting surface

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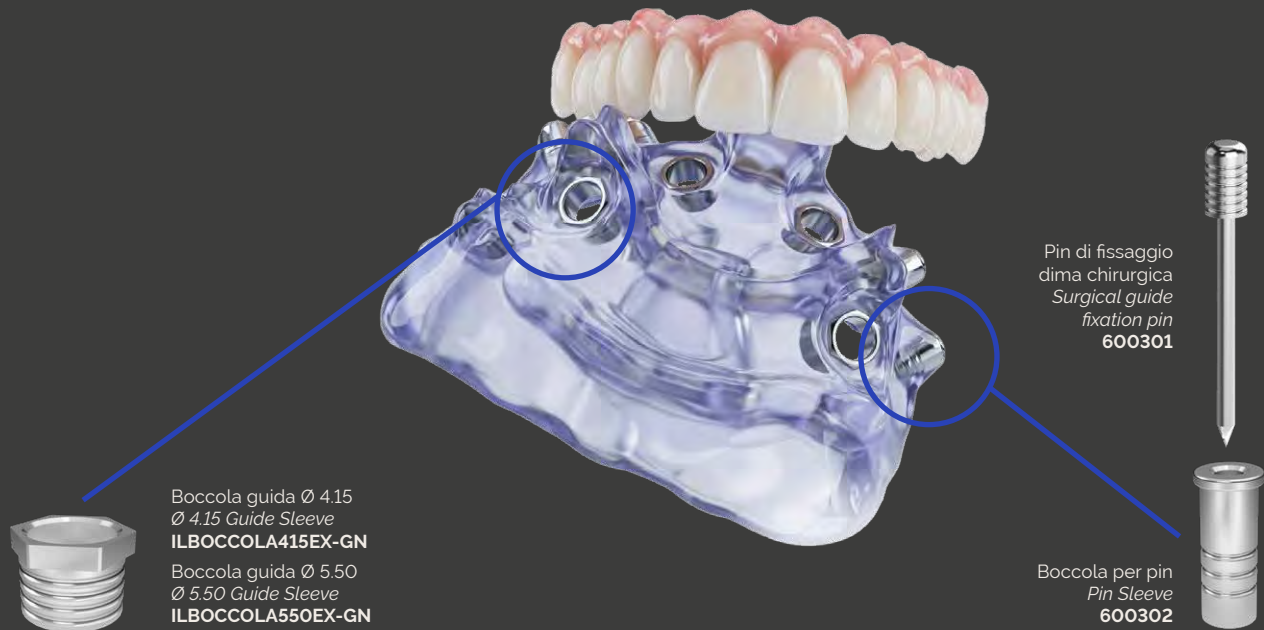
GUIDE SLEEVES

The guide sleeves are hexagonal-headed cylinders (indexed) that are incorporated into the surgical guide. Their purpose is to guide the instruments during surgery. By ensuring a fixed ratio of 9.00 mm between the instrument contact point on the bushing and the implant portion to be inserted into the bone, the planning created with the 3P GUIDE software is respected.

The guide sleeves are 4.00 mm high and come in two different diameters. The LARGE sleeve has an internal diameter of 5.50 mm, while the SMALL sleeve has an internal diameter of 4.15 mm. The SMALL sleeve is only suitable for implants with a diameter of 3.40 mm or 3.80 mm.

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PIN SLEEVES

Pin sleeves are metal cylinders incorporated into the surgical guide and are used to stabilize the guide, especially in cases where the guide cannot support sufficient residual elements. The guide must be stable to ensure safety during surgery. The absence of micromovements allows the desired accuracy to be achieved.

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SEQUENZA PASSAGGI / STEP SEQUENCE	H 6.0	H 8.5	H 10.0	H 11.5	H 13.0	H 15.0
	1. tissue punch	1. tissue punch	1. tissue punch	1. tissue punch	1. tissue punch	1. tissue punch
2.	crestal mill 3,40	crestal mill 3,40	crestal mill 3,40	crestal mill 3,40	crestal mill 3,40	crestal mill 3,40
3.	initial drill	initial drill	initial drill	initial drill	initial drill	initial drill
SEQUENZA FRESE / DRILLS SEQUENCE (Ø x H)	2x6	2x6	2x6	2x6	2x6	2x6
	3,4x6	2x8,5	2x8,5	2x8,5	2x8,5	2x8,5
		3,4x8,5	2x10	2x10	2x10	2x10
			3,4x10	2x11,5	2x11,5	2x11,5
				3,4x10	2x13	2x13
				3,4x11,5	3,4x10	2x15
					3,4x11,5	3,4x10
					3,4x13	3,4x11,5
						3,4x13
						3,4x15

Ø 3.80 mm 

MOUNTER COLLO / NECK IL® 1,8 - 2,8 - 3,8 / IL WIDE® 0,9 - 1,8 - 2,8
TIPO DI BOCCOLA / SLEEVE TYPE 4,15 - 5,50

SEQUENZA PASSAGGI / STEP SEQUENCE	H 6.0	H 8.5	H 10.0	H 11.5	H 13.0	H 15.0
	1. tissue punch	1. tissue punch	1. tissue punch	1. tissue punch	1. tissue punch	1. tissue punch
2.	crestal mill 3,80	crestal mill 3,80	crestal mill 3,80	crestal mill 3,80	crestal mill 3,80	crestal mill 3,80
3.	initial drill	initial drill	initial drill	initial drill	initial drill	initial drill
SEQUENZA FRESE / DRILLS SEQUENCE (Ø x H)	2x6	2x6	2x6	2x6	2x6	2x6
	3,4x6	2x8,5	2x8,5	2x8,5	2x8,5	2x8,5
	3,8x6	3,4x8,5	2x10	2x10	2x10	2x10
		3,8x8,5	3,4x10	2x11,5	2x11,5	2x11,5
			3,8x10	3,4x10	2x13	2x13
				3,4x11,5	3,4x10	2x15
				3,8x10	3,4x11,5	3,4x10
				3,8x11,5	3,4x13	3,4x11,5
					3,8x10	3,4x13
					3,8x11,5	3,4x15

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The protocol must be customized according to the bone characteristics and the type of surgical procedure to be performed.

Very dense bone (D1-D2): it is preferable to use all drills to avoid bone stress.

Soft bone (D3-D4): consider under-preparing the site in diameter to achieve greater primary stability.

Ø 4.25 mm 

MOUNTER COLLO / NECK IL® 18 - 2,8 - 3,8 / IL WIDE® 0,9 - 1,8 - 2,8
TIPO DI BOCCOLA / SLEEVE TYPE 5,50

SEQUENZA PASSAGGI / STEP SEQUENCE	H 6.0 H 8.5 H 10.0 H 11.5 H 13.0 H 15.0					
	1.	2.	3.			
	tissue punch	tissue punch	tissue punch	tissue punch	tissue punch	tissue punch
	crestal mill 4,25	crestal mill 4,25	crestal mill 4,25	crestal mill 4,25	crestal mill 4,25	crestal mill 4,25
	initial drill	initial drill	initial drill	initial drill	initial drill	initial drill
SEQUENZA FRESE / DRILLS SEQUENCE (Ø x H)	2x6	2x6	2x6	2x6	2x6	2x6
	3,4x6	2x8,5	2x8,5	2x8,5	2x8,5	2x8,5
	3,8x6	3,4x8,5	2x10	2x10	2x10	2x10
	4,25x6	3,8x8,5	3,4x10	2x11,5	2x11,5	2x11,5
		4,25x8,5	3,8x10	3,4x10	2x13	2x13
			4,25x10	3,4x11,5	3,4x10	2x15
				3,8x10	3,4x11,5	3,4x10
				3,8x11,5	3,4x13	3,4x11,5
				4,25x10	3,8x10	3,4x13
				4,25x11,5	3,8x11,5	3,4x15
				3,8x13	3,8x10	
				4,25x10	3,8x11,5	
				4,25x11,5	3,8x13	
				4,25x13	3,8x15	
					4,25x10	
					4,25x11,5	
					4,25x13	
					4,25x15	

Ø 5.00 mm 

MOUNTER COLLO / NECK IL® 2,8 - 3,8 / IL WIDE® 0,9 - 1,8 - 2,8
TIPO DI BOCCOLA / SLEEVE TYPE 5,50

SEQUENZA PASSAGGI / STEP SEQUENCE	H 6.0 H 8.5 H 10.0 H 11.5 H 13.0 H 15.0					
	1.	2.	3.			
	tissue punch	tissue punch	tissue punch	tissue punch	tissue punch	tissue punch
	crestal mill 5,00	crestal mill 5,00	crestal mill 5,00	crestal mill 5,00	crestal mill 5,00	crestal mill 5,00
	initial drill	initial drill	initial drill	initial drill	initial drill	initial drill
SEQUENZA FRESE / DRILLS SEQUENCE (Ø x H)	2x6	2x6	2x6	2x6	2x6	2x6
	3,4x6	2x8,5	2x8,5	2x8,5	2x8,5	2x8,5
	3,8x6	3,4x8,5	2x10	2x10	2x10	2x10
	4,25x6	3,8x8,5	3,4x10	2x11,5	2x11,5	2x11,5
	5x6	4,25x8,5	3,8x10	3,4x10	2x13	2x13
		5x8,5	4,25x10	3,4x11,5	3,4x10	2x15
			5x10	3,8x10	3,4x11,5	3,4x10
				3,8x11,5	3,4x13	3,4x11,5
				4,25x10	3,8x10	3,4x13
				4,25x11,5	3,8x11,5	3,4x15
				5x10	3,8x13	3,8x10
				5x11,5	4,25x10	3,8x11,5
					4,25x11,5	3,8x13
					4,25x13	3,8x15
					5x10	4,25x10
					5x11,5	4,25x11,5
					5x13	4,25x13
						4,25x15
						5x10
						5x11,5
						5x13
						5x15

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Quick&Easy is a prosthetic protocol based on the intramucosal positioning of the IL® and IL WIDE® implants to create single crowns, bridges, or full arches, transferring the B.O.P.T.* concept to the straight abutment: Iconic Abutment (IA-BOPT®).

The **Quick&Easy** procedure is based on the use of a few components contained in the Iconic Box and, for some types of full-arch restorations, on the milled titanium bar that is screwed into the dental prosthetic elements.

Biography

Drs. Giuseppe and Costantino are the owners of the Vignato Clinic in Vicenza and they began developing a clinical protocol for convergent-neck implants in 2013. Since 2025, they have been part of the Scientific Committee of Leader Medica, contributing to the company's scientific and cultural development.

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The Quick & Easy brand is owned by Dr. Costantino Vignato, with exclusive use by Leader Medica Srl.

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ICONIC ENGAGING HEX ABUTMENT
FOR SINGLE CROWNS

H6	ILTB6N	H6	ILWTB6N
H8	ILTB8N		
H10	ILTB10N		

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ICONIC ROTATING ABUTMENT
FOR BRIDGES

H6	ILTB6R	H6	ILWTB6R
H8	ILTB8R		
H10	ILTB10R		

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H2	ILVGDR2	H2	ILWVGDR2
H3	ILVGDR3	H3	ILWVGDR3
H5	ILVGDR5	H5	ILWVGDR5
H7	ILVGDR7		

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TEMPORARY ABUTMENTS FOR
PASSIVATION TESTS USED IN
FULL-ARCH RESTORATIONS

ILMPN	ILWMPN
ILMPR	ILWMPR
	ILWMPAN
	ILWMPAR

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CYLINDER HEAD PROSTHETIC
SCREW FOR FLAT SEAL
ILVPR

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SHORT AND LONG MANUAL DRIVER
HEXAGON 1.28

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TORX CYLINDER HEAD PROSTHETIC
SCREW FOR FLAT SEAL
ILVPRTX

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The Iconic Abutment (IA-BOPT[®]) replicates the anatomical shape of a natural tooth and allows for the creation of any type of prosthesis with a single component. The vertical preparation allows the crown to be shaped so that it perfectly adapts to the gingival profile, improving the aesthetic and functional integration in the mouth.

When the implant is off-axis, the access hole for the screw may be in an unaesthetic position. For this reason, the Iconic Box includes a dynamic screw, which allows for angled screwing up to 30°, ensuring the hole exits in the ideal position, optimizing aesthetics and practicality.

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*“The prosthesis that leaves
room for soft tissue”*

(A-B.O.P.T.®

Protocollo Protesico Iconic Abutment

Iconic Abutment Prosthetic Protocol

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PLACEMENT OF THE INTRAMUCOSAL IMPLANT

The convergent implant facilitates the design of the emerging prosthesis, allowing for greater space for soft tissue and a better aesthetic result.

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HEALING SCREW

Insert the healing screw immediately after implant placement. The straight healing screw is the smallest crown that can be fabricated.

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Iconic rotating abutment

FULL ARCH:

Iconic Abutment (rotating) + Milled bar + components for intermediate try-ins

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IMPRESSION TAKING

Once the tissues have healed, generally after 2-3 months, depending on the condition of the implant site.

- Customized transfers (if CAD/CAM) or standard (if analog)
- Intraoral scanning (preferred for precision and biological savings)

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CROWN DESIGN ON THE ICONIC ABUTMENT

Digital CAD design with modeling of the crown directly on the abutment.

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SPLASH EFFECT

When we remove the healing screw and insert the cemented crown onto the iconic abutment, we design the shape of the soft tissue.

Si consiglia coppia di serraggio protesica a 25 N·cm | *The advised prosthetic tightening torque is 25 N·cm*

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HEALING SCREWS

Designed to ensure guided and harmonious gingival shaping, healing screws come in different heights and diameters to adapt to the individual mucosal profile. The cylindrical-convergent shape favors the thickness of the soft tissue and optimally prepares the site for the future prosthetic component.

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ANALOGS AND TRANSFERS

Implant analogs and analog impression transfers (pick-up and tear-off) allow for an accurate replication of the implant's three-dimensional position within the model. Their dimensional accuracy is crucial to ensure accuracy during the laboratory phase. Leader Medica recommends always using disposable transfers and analogs to ensure maximum precision in the implant connection.

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SCANBODY AND DIGITAL ANALOGS

Optical impression scanbodies and digital analogs accurately reproduce implant position and angle in CAD/CAM workflows. They are manufactured with the highest precision to ensure the highest fidelity in transmitting clinical data.

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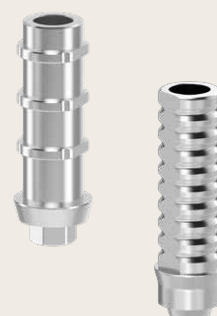
TEMPORARY ABUTMENTS

Temporary abutments provide for the optimal management of the emergence profile during the healing and tissue conditioning phase. They allow for immediate or delayed loading techniques with stability, simplicity, and predictability. They are available in a non-rotating version, suitable for single-implant restorations, and a rotating version, recommended for multiple-implant prosthetic rehabilitations, such as bridges.

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STRAIGHT AND ANGLED ABUTMENTS

High-precision medical-grade titanium prosthetic solutions, designed for cemented prosthetics in single or multiple restorations, in both front and posterior areas. Being straight or angled, the abutments ensure passive fit, strength, and reliability.

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CASTABLE AND COBALT-CHROME CANNULAE

Castable and CoCr cannulae are ideal for creating customized prosthetic structures by casting or milling on screw-retained prostheses. Strong and moldable, they ensure excellent precision in the creation of custom-made superstructures. Compatible with CAD-CAM and traditional techniques.

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LOCATOR

The use of Locators with metal caps allows for a secure and versatile connection of removable dentures, providing a stable attachment system for optimal comfort and greater chewing efficiency.

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MUA

MUA (multi-unit abutments), designed for all-on-4/6 implants, allow for angulation correction, ensuring a functional and well-fitting prosthesis.

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HEALING SCREWS WITH CORONAL ENLARGEMENT

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PER IMPRONTA ANALOGICA FOR ANALOGIC PRINT

ILTR

Transfer pick up emergenza dritta IL®
Transfer pick up straight emergence IL®

ILAN3828

Analogo 3828 IL®
Analogue 3828 IL®

ILTRPULL

Transfer pull up in peek IL®
Transfer pull up in peek IL®

Transfer ILTR: vite ILVTR inclusa / screw ILVTR included

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Analogo digitale / Digital analogue ILANCADxxyy: vite ILANSCREW inclusa / screw ILANSCREW included

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Moncone / Abutment ILMD h6 - h8: vite ILVPRSV inclusa / screw ILVPRSV included
Moncone / Abutment ILMD h10 - h15: vite ILVPR inclusa / screw ILVPR included

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ILTBF6N

Iconic Abutment h6 non rotante, B.O.P.T.®, IL®
Iconic Abutment h6, B.O.P.T.®, non-rotating, IL®

ILTBF6R

Iconic Abutment h6 rotante, B.O.P.T.®, IL®
Iconic Abutment h6, B.O.P.T.®, rotating, IL®

ILTBF8N

Iconic Abutment h8 non rotante, B.O.P.T.®, IL®
Iconic Abutment h8, B.O.P.T.®, non-rotating, IL®

ILTBF8R

Iconic Abutment h8 rotante, B.O.P.T.®, IL®
Iconic Abutment h8, B.O.P.T.®, rotating, IL®

ILTBF10N

Iconic Abutment h10 non rotante, B.O.P.T.®, IL®
Iconic Abutment h10, B.O.P.T.®, non-rotating, IL®

ILTBF10R

Iconic Abutment h10 rotante, B.O.P.T.®, IL®
Iconic Abutment h10, B.O.P.T.®, rotating, IL®

Iconic Abutment ILTBFxx: vite ILVPR inclusa / screw ILVPR included

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T-Base finestrata / windowed ILTBFFxx: vite ILVPR inclusa / screw ILVPR included

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T-Base Biologic ILTBBLxx: vite ILVPRSV inclusa / screw ILVPRSV included

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T-Base Biologic Finestrata ILTBFLxx: vite ILVPR inclusa / screw ILVPR included

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T-Base with 1.5 intramucosal height and h4 cementing cone, rotating, IL®

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T-Base with 1.5 intramucosal height and h6 cementing cone, rotating, IL®

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T-Base with 2.5 intramucosal height and h4 cementing cone, rotating, IL®

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T-Base with 2.5 intramucosal height and h6 cementing cone, rotating, IL®

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T-Base with 3.5 intramucosal height and h4 cementing cone, rotating, IL®

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T-Base with 3.5 intramucosal height and h6 cementing cone, rotating, IL®

T-Base ILTBxx: vite ILVPR inclusa / screw ILVPR included

SU RICHIESTA / ON REQUEST: vite ILVPRTX / screw ILVPRTX

BASETTE SISTEMI CAD/CAM DIGITALI DIGITAL CAD/CAM SYSTEMS BASES

ILTBCER05

T-Base h0,5 compatibile con sistemi CAD/CAM digitali Cerec®, IL®

T-Base h0.5 compatible with digital Cerec® CAD/CAM systems IL®

ILTBCER15

T-Base h1,5 compatibile con sistemi CAD/CAM digitali Cerec®, IL®

T-Base h1.5 compatible with digital Cerec® CAD/CAM systems IL®

T-Base Cer ILTBFCERxx: vite ILVPR inclusa / screw ILVPR included

ILTBBL CER15

T-Base Biologic h1,5 compatibile con sistemi CAD/CAM digitali Cerec®, IL®

Biologic T-Base h1.5 compatible with digital Cerec® CAD/CAM systems IL®

T-Base Biologic Cer ILTBBL CER15: vite ILVPRSV inclusa / screw ILVPRSV included

T-BASE IN CROMO COBALTO CON CANNULA CALCINABILE PER PROTESI AVVITATA COBALT-CHROME T-BASE WITH CALCINABLE CANNULA FOR SCREW-RETAINED PROSTHESIS

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Cobalt-Chrome T-Base with calcinable cannula, non-rotating, IL®

ILTBCRCOR

T-Base rotante in Cromo Cobalto con cannula calcinabile IL®

Cobalt-Chrome T-Base with calcinable cannula, rotating, IL®

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MUA - TRANSFER AND ANALOGUES

IL900003

Transfer pick up per MUA
Transfer pick up for MUA

TIS900001

Analogo per MUA
Analogue for MUA

Transfer pick up IL900003: vite IL900004 inclusa / screw IL900004 included

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Scanbody ILINTRA: vite IL100003 inclusa / screw IL100003 included

Analogo CAD /CAD analogue ILANCADMUA: vite TISANSCREW128 inclusa / screw TISANSCREW128 included

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Moncone provvisorio /Temporary abutment IL100222: vite IL100003 inclusa / screw IL100003 included

MUA - MONCONE CALCINABILE

MUA - CALCINABLE ABUTMENT

IL900002

Moncone calcinabile per MUA
Calcinable abutment for MUA

Moncone calcinabile /Calcinable abutment IL900002: vite IL100003 inclusa / screw IL100003 included

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IL900004

Vite transfer per MUA
Transfer screw for MUA

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The advised prosthetic tightening torque is 25 N·cm

IL100003:

The advised prosthetic tightening torque is 15 N·cm

Strumentario IL[®]

IL[®], Instruments

STRUMENTARIO PER COMPONENTI PROTESICHE INSTRUMENTS FOR PROSTHETIC COMPONENTS

TCOTORXC

Cacciavite Torx Corto
Short Torx screwdriver

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The convergent evolution takes shape through the combination of IL WIDE® implant and the UNICO® abutment, specifically designed to generate the intramucosal convergent effect on this implant line, but applicable across all so-called bone level Leader Medica implant lines, regardless of the connection type.

UNICO® is designed with an interface portion to ensure stable fit with the fixture, a intramucosal portion shaped to recreate an intramucosal emerging profile, and an upper portion suitable for connection to a prosthetic superstructure.

Its geometry, designed to preserve the biological width according to the principles of B.O.P.T.®, allows bone level osseointegrated implants to be transformed into intramucosal prosthetic configurations without having to surgically replace the implant. This results in a perfect solution for promoting soft tissue healing, stabilizing the mucosal seal, and optimizing the distribution of masticatory forces.

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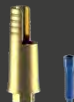
Unico Abutment straight h6, non-rotating, B.O.P.T.®, IL WIDE®

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Unico Abutment straight h8, non-rotating, B.O.P.T.®, IL WIDE®

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Unico Abutment straight h10, non-rotating, B.O.P.T.®, IL WIDE®

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Unico Abutment 15° angled h6, non-rotating, B.O.P.T.®, IL WIDE®

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Unico Abutment straight h6, rotating, B.O.P.T.®, IL WIDE®

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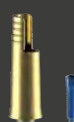
Unico Abutment straight h8, rotating, B.O.P.T.®, IL WIDE®

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Unico Abutment straight h10, rotating, B.O.P.T.®, IL WIDE®

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Unico Abutment 15° angled h6, rotating, B.O.P.T.®, IL WIDE®

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*One prosthetic component for
each implant connection
from IL WIDE® to MITO, FIVE...*

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/ FREEDOM OF DESIGN

The UNICO® abutment follows the design principles of Dr. Ignazio Loi's B.O.P.T.® preparation, giving the prosthodontist complete freedom to determine the vertical position of the prosthetic margin closure — also in relation to the depth at which the implant is placed.

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/ FOR EVERY HEIGHT

The UNICO® system is modular and allows the clinician to choose the most suitable intramucosal height according to the patient's gingival biotype.

It is compatible with both screw-retained and cement-retained prostheses, offering a standardized abutment available in different heights to adapt to various clinical situations.

BENEFITS

- Compatible with both **digital and analog workflows**
- Recommended for cases with **severely angled or non-parallel implants**
- Available **in straight or angled version**, suitable for the prosthetic restoration of all bone level implants
- Individualized **esthetic closure** on the cylindrical portion

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SPECIFICATIONS

- Cementation taper: **6°**
- Angulation: **15°**
- Versions: **rotating and non-rotating**, for single restorations or multi-unit bridges

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TRANSFER, ANALOGHI E SCANBODY TRANSFER, ANALOGUES AND SCANBODY

PER IMPRONTA ANALOGICA FOR ANALOGIC PRINT

ILWTR

Transfer pick up emergenza dritta IL WIDE[®]
Transfer pick up straight emergency IL WIDE[®]

ILWAN3809

Analogo 3809 IL WIDE[®]
Analogue 3809 IL WIDE[®]

ILWTRST

Transfer a strappo emergenza anatomica IL WIDE[®]
Closed tray transfer anatomical emergency IL WIDE[®]

Transfer ILWTR - ILWTRST: vite ILVTR inclusa / screw ILVTR included

PER IMPRONTA DIGITALE FOR DIGITAL PRINT

ILWIOS

Scanbody IL WIDE®
Scanbody IL WIDE®

ILWANCAD4218

Analogo digitale Ø 4,25 Collo 1,8 IL WIDE®
Digital analogue Ø 4.25 Neck 1.8 IL WIDE®

ILWIOS-2

Scanbody IL WIDE® Slim
Scanbody IL WIDE® Slim

ILWANCAD4228

Analogo digitale Ø 4,25 Collo 2,8 IL WIDE®
Digital analogue Ø 4.25 Neck 2.8 IL WIDE®

ILWANCAD3809

Analogo digitale Ø 3,80 Collo 0,9 IL WIDE®
Digital analogue Ø 3.80 Neck 0.9 IL WIDE®

ILWANCAD5009

Analogo digitale Ø 5,00 Collo 0,9 IL WIDE®
Digital analogue Ø 5.00 Neck 0.9 IL WIDE®

ILWANCAD3818

Analogo digitale Ø 3,80 Collo 1,8 IL WIDE®
Digital analogue Ø 3.80 Neck 1.8 IL WIDE®

ILWANCAD5018

Analogo digitale Ø 5,00 Collo 1,8 IL WIDE®
Digital analogue Ø 5.00 Neck 1.8 IL WIDE®

ILWANCAD3828

Analogo digitale Ø 3,80 Collo 2,8 IL WIDE®
Digital analogue Ø 3.80 Neck 2.8 IL WIDE®

ILWANCAD5028

Analogo digitale Ø 5,00 Collo 2,8 IL WIDE®
Digital analogue Ø 5.00 Neck 2.8 IL WIDE®

ILWANCAD4209

Analogo digitale Ø 4,25 Collo 0,9 IL WIDE®
Digital analogue Ø 4.25 Neck 0.9 IL WIDE®

ILWANCAD6028

Analogo digitale Ø 6,00 Collo 2,8 IL WIDE®
Digital analogue Ø 6.00 Neck 2.8 IL WIDE®

Scanbody: vite ILVPR inclusa / screw ILVPR included

Analogo digitale / Digital analogue ILWANCADxxyy: vite ILANSCREW inclusa / screw ILANSCREW included

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Pilastri provvisori / Temporary abutments: vite ILVPR inclusa / screw ILVPR included

MONCONI DRITTI PER PROTESI CEMENTATA/AVVITATA STRAIGHT ABUTMENTS FOR CEMENTED/SCREW-RETAINED PROSTHESES

		ILWMD6 Moncone dritto Profile non rotante h6 IL WIDE® <i>Profile straight abutment non-rotating h6 IL WIDE®</i>			ILWMD6R Moncone dritto Profile rotante h6 IL WIDE® <i>Profile straight abutment rotating h6 IL WIDE®</i>
		ILWMD8 Moncone dritto Profile non rotante h8 IL WIDE® <i>Profile straight abutment non-rotating h8 IL WIDE®</i>			ILWMD8R Moncone dritto Profile rotante h8 IL WIDE® <i>Profile straight abutment rotating h8 IL WIDE®</i>
		ILWMD10 Moncone dritto Profile non rotante h10 IL WIDE® <i>Profile straight abutment non-rotating h10 IL WIDE®</i>			ILWMD10R Moncone dritto Profile rotante h10 IL WIDE® <i>Profile straight abutment rotating h10 IL WIDE®</i>
		ILWMD15 Moncone dritto Profile non rotante h15 IL WIDE® <i>Profile straight abutment non-rotating h15 IL WIDE®</i>			

Moncone / Abutment ILWMD h6 - h8: vite ILVPRSV inclusa / screw ILVPRSV included
Moncone / Abutment ILWMD h10 - h15: vite ILVPR inclusa / screw ILVPR included

MONCONI ANGOLATI PER PROTESI CEMENTATA ANGLED ABUTMENTS FOR CEMENTED PROSTHESES

		ILWMA15°1 Moncone angolato Profile 15° h1 IL WIDE® <i>Profile 15° angled abutment h1 IL WIDE®</i>			ILWMA25°1 Moncone angolato Profile 25° h1 IL WIDE® <i>Profile 25° angled abutment h1 IL WIDE®</i>
		ILWMA15°2 Moncone angolato Profile 15° h2 IL WIDE® <i>Profile 15° angled abutment h2 IL WIDE®</i>			ILWMA25°2 Moncone angolato Profile 25° h2 IL WIDE® <i>Profile 25° angled abutment h2 IL WIDE®</i>

Moncone / Abutment ILWMAxx°y: vite ILVPR inclusa / screw ILVPR included

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ICONIC ABUTMENT FINESTRATO PER B.O.P.T.® WINDOWED ICONIC ABUTMENT FOR B.O.P.T.®

		ILWTBF6N Iconic Abutment h6 non rotante, B.O.P.T.®, IL WIDE® <i>Iconic Abutment h6, B.O.P.T.®, non-rotating, IL WIDE®</i>			ILWTBF6R Iconic Abutment h6 rotante, B.O.P.T.®, IL WIDE® <i>Iconic Abutment h6, B.O.P.T.®, rotating, IL WIDE®</i>
		ILWTBF8N Iconic Abutment h8 non rotante, B.O.P.T.®, IL WIDE® <i>Iconic Abutment h8, B.O.P.T.®, non-rotating, IL WIDE®</i>			ILWTBF8R Iconic Abutment h8 rotante, B.O.P.T.®, IL WIDE® <i>Iconic Abutment h8, B.O.P.T.®, rotating, IL WIDE®</i>
		ILWTBF10N Iconic Abutment h10 non rotante, B.O.P.T.®, IL WIDE® <i>Iconic Abutment h10, B.O.P.T.®, non-rotating, IL WIDE®</i>			ILWTBF10R Iconic Abutment h10 rotante, B.O.P.T.®, IL WIDE® <i>Iconic Abutment h10, B.O.P.T.®, rotating, IL WIDE®</i>

Iconic Abutment ILWTBFxx: vite ILVPR inclusa / screw ILVPR included

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T-Base Dyn windowed with platform-level closure and h4 cementing cone, non-rotating, IL WIDE®

T-Base finestrata / windowed ILWTBFxx: vite ILVPR inclusa / screw ILVPR included

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T-Base Dyn windowed with platform-level closure and h4 cementing cone, rotating, IL WIDE®

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T-Base with platform-level closure and h4 cementing cone, non-rotating, IL WIDE®

T-Base ILWTBxx: vite ILVPR inclusa / screw ILVPR included

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T-Base with platform-level closure and h6 cementing cone, non-rotating, IL WIDE®

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T-Base with platform-level closure and h4 cementing cone, rotating, IL WIDE®

T-Base ILWTBxx: vite ILVPR inclusa / screw ILVPR included

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T-Base with platform-level closure and h6 cementing cone, rotating, IL WIDE®

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NON-ROTATING BASES WITH INTRAMUCOSAL H

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T-Base ILWTBxx: vite ILVPR inclusa / screw ILVPR included

SU RICHIESTA / ON REQUEST: vite ILVPRTX / screw ILVPRTX

BASETTE ROTANTI CON H INTRAMUCOSA

ROTATING BASES WITH INTRAMUCOSAL H

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T-Base ILWTBxx: vite ILVPR inclusa / screw ILVPR included

SU RICHIESTA / ON REQUEST: vite ILVPRTX / screw ILVPRTX

BASETTE SISTEMI CAD/CAM DIGITALI

DIGITAL CAD/CAM SYSTEMS BASES

ILWTBCER15

T-Base h1,5 compatibile con sistemi CAD/CAM digitali Cerec®, IL WIDE®
T-Base h1.5 compatible with Cerec® digital CAD/CAM systems IL WIDE®

T-Base Cer ILTBCERxx: vite ILVPR inclusa / screw ILVPR included

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COBALT-CHROME T-BASE WITH CALCINABLE CANNULA FOR SCREW-RETAINED PROSTHESIS

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Cobalt-Chrome T-Base with calcinable cannula, non-rotating, IL WIDE®

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Cobalt-Chrome T-Base with calcinable cannula, rotating, IL WIDE®

T-Base ILWTBCRCON, ILWTBCRCOR: vite ILVPR inclusa / screw ILVPR included

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MUA angolato / Angled MUA: vite ILVPRMUA inclusa / screw ILVPRMUA included

MUA T-BASE CROMO COBALTO CON CANNULA CALCINABILE

COBALT-CHROME T-BASE MUA WITH CALCINABLE CANNULA

ILBACCR
T-Base in Cromo Cobalto con cannula
calcinabile per MUA
*Cobalt-Chrome T-Base with calcinable cannula
for MUA*

T-Base MUA: vite IL100003 inclusa / screw IL100003 included

MUA T-BASE TITANIO

TITANIUM T-BASE MUA

ILBATI
T-Base in titanio per MUA
Titanium T-Base for MUA

ILBATISH
T-Base short in titanio per MUA
Titanium T-Base short for MUA

T-Base MUA: vite IL100003 inclusa / screw IL100003 included

SU RICHIESTA / ON REQUEST:
vite TIS100003TX / screw TIS100003TX

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The advised prosthetic tightening torque is 25 N-cm

IL100003:

The advised prosthetic tightening torque is 15 N-cm

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IL[®] Zygomatic, Nasal and Pterygoid Implant

IL[®] ZIGOMATICO Ø 3,75 IL[®] ZYGOMATIC Ø 3.75

h 35 mm	ILZY3735	ILZY4235
h 37,5 mm	ILZY3737	ILZY4237
h 40 mm	ILZY3740	ILZY4240
h 42,5 mm	ILZY3742	ILZY4242
h 45 mm	ILZY3745	ILZY4245
h 47,5 mm	ILZY3747	ILZY4247
h 50 mm	ILZY3750	ILZY4250
h 52,5 mm	ILZY3752	ILZY4252
h 55 mm	ILZY3755	ILZY4255
h 60 mm	ILZY3760	ILZY4260

IL[®] NASALE Ø 4,00 IL[®] NASAL Ø 4.00

h 24 mm	ILNS4024
h 26 mm	ILNS4026
h 28 mm	ILNS4028
h 30 mm	ILNS4030

IL[®] PTERIGOIDEO Ø 4,00 IL[®] PTERYGOID Ø 4.00

h 18 mm	ILPT4018
h 20 mm	ILPT4020
h 22 mm	ILPT4022

IL[®] Zigomatico, Nasale e Pterigoideo

Componenti Protetiche

IL[®] Zygomatic, Nasal and Pterygoid, Prosthetic Components

MUA – MONCONI MULTI UNIT - MULTI UNIT ABUTMENT

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The indications and images are just general instructions that should be adapted on a case-by-case basis according to the patient's needs. These recommendations are not intended to replace specialized training or the expertise of the clinician, who—also based on personal experience—can make different decisions.

ESSENTIAL PROCEDURES TO ENSURE THE SAFETY, EFFECTIVENESS, AND FUNCTIONALITY OF INSTRUMENTS

All surgical instruments must be thoroughly cleaned, disinfected, and sterilized before first use and before each subsequent reuse. It is essential to check the proper functioning of the instruments before each procedure and, if necessary, to immediately replace any instruments showing signs of wear or malfunction. Failure to follow these instructions may increase the risk of infection and intraoperative complications for the patient.

PREPARING THE AREA FOR IMPLANT PLACEMENT

Every procedure must be performed gradually and carefully, avoiding sudden movements that could compromise safety. For traditional surgery check the correct positioning of the drill stops, ensuring they are fully tightened. Choose the length in accordance with the implant insertion plan and the given warnings. Periodically check surgical micromotors used with drills and taps to prevent malfunctions. Regularly check the drills for wear and replace them if they lose effectiveness in order to ensure precision and safety during implant site preparation.

PROSTHESIS MAINTENANCE

To prevent potential mechanical and biological complications associated with prosthetic implants—which are widely documented in the literature—it is essential for the patient to maintain proper oral hygiene and undergo scheduled periodic checkups. These measures help prolong the functional life of the device. Periodic checks of the calibration of abutment and prosthetic screws are also recommended to ensure the stability and efficiency of the device. If a feeling of instability or malfunction occurs, the patient should be advised to seek prompt medical attention.

DISPOSING OF INSTRUMENTS AND IMPLANTS

Surgical instruments, consisting of small metal parts, can be disposed of as metal waste only if properly cleaned. Otherwise, they must be treated as biological waste. Similarly, dental implants removed from the patient's mouth must be treated as biological waste. For each disposal procedure, it is recommended to refer to the local regulations in force regarding the management of medical waste.

BIBLIOGRAFIA

Bibliography

FMC (FIBROBLAST MICRO CHANNEL)

Harrison, Ross Granville. "The cultivation of tissues in extraneous media as a method of morphogenetic study". *The Anatomical Record*. (1912) 6 (4): 181-193

Weiss, Paul. "Experiments on cell and axon orientation in vitro: The role of colloidal exudates in tissue organization". *Journal of Experimental Zoology*. (1945) 100 (3): 353-386.

Brunette DM. The effects of implant surface topography on the behavior of cells. *Int J Oral Maxillofac Implants*. 1988 Winter;3(4):231-46. PMID: 3075965.

Ricci JL, Grew JC, Alexander H. Connective-tissue responses to defined biomaterial surfaces. I. Growth of rat fibroblast and bone marrow cell colonies on microgrooved substrates. *J Biomed Mater Res A*. 2008 May;85(2):313-25. doi: 10.1002/jbma.31379. PMID: 17688290.

Nevins M, Nevins ML, Camelo M, Boyesen JL, Kim DM. Human histologic evidence of a connective tissue attachment to a dental implant. *Int J Periodontics Restorative Dent*. 2008 Apr;28(2):111-21. PMID: 18546807.

Bria J, Guzzo M, Morelli S, De Santo MP, Chiarella E, Filardo A, Barone S, De Bartolo L, Covani U, Giudice A, Di Vito A. - In Vitro Cellular Responses to Titanium Surface Topographies: Toward Improved Mucointegration of Dental Implants - Submission in progress.

HPA (HIGH PERFORMANCE ANODIZING)

Corvino E, Pesce P, Mura R, Marcano E, Canullo L. Influence of Modified Titanium Abutment Surface on Peri-implant Soft Tissue Behavior: A Systematic Review of In Vitro Studies. *Int J Oral Maxillofac Implants*. 2020 May/Jun;35(3):503-519. doi: 10.11607/jomi.8110. PMID: 32406646.

Traver-Méndez V, Camps-Font O, Ventura F, Nicolau-Sansó MA, Subirà-Pifarré C, Figueiredo R, Valmaseda-Castellón E. In Vitro Characterization of an Anodized Surface of a Dental Implant Collar and Dental Abutment on Peri-Implant Cellular Response. *Materials* (Basel). 2023 Sep 1;16(17):6012. doi: 10.3390/ma16176012. PMID: 37687706; PMCID: PMC10489139.

Tingting Wang, Lina Wang, Qianqian Lu, Zhen Fan, Influence of anodized titanium abutments on the esthetics of the peri-implant soft tissue: A clinical study. *The Journal of Prosthetic Dentistry*, Volume 125, Issue 3, 2021, Pages 445-452, ISSN 0022-3913, <https://doi.org/10.1016/j.prosdent.2019.12.023>.

Farrag KM, Khamis MM. Effect of anodized titanium abutment collars on peri-implant soft tissue: A split-mouth clinical study. *J Prosthet Dent*. 2023 Jul;130(1):59-67. doi: 10.1016/j.prosdent.2021.09.019. Epub 2021 Oct 29. PMID: 34756606.

Liu W, Gregory RL, Yang CC, Hamada Y, Lin WS. The effects of anodization and instrumentation on titanium abutment surface characteristics and biofilm formation. *J Prosthodont*. 2025 Jan 15. doi: 10.1111/jopr.14009. Epub ahead of print. PMID: 39810607.

Stefano Gracis, Nitzan Bichacho, Arturo Llobell, Leila Jahangiri. The Influence of Implant Neck Features and Abutment Diameter on Hard and Soft Tissues Around Single Implants Placed in Healed Ridges: Clinical Criteria for Selection. Article in *The International journal of periodontics & restorative dentistry*, January 2020. DOI: 10.11607/prd.4151.

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DAS (DOUBLE ACIDIFIED SURFACE)

Jang TS, Jung HD, Kim S, Moon BS, Baek J, Park C, Song J, Kim HE. Multiscale porous titanium surfaces via a two-step etching process for improved mechanical and biological performance. *Biomed Mater*. 2017 Mar 15;12(2):025008. doi: 10.1088/1748-605X/aa5d74. PMID: 28296644.

Giner L, Mercadé M, Torrent S, Punset M, Pérez RA, Delgado LM, Gil FJ. Double acid etching treatment of dental implants for enhanced biological properties. *J Appl Biomater Funct Mater*. 2018 Apr;16(2):83-89. doi: 10.5301/jabfm.5000376. PMID: 28885666.

Santos Marino J, Cortés-Bretón Brinkmann J, García-Gil I, Martínez-Rodríguez N, Fraile JF, Barona Dorado C, Martínez-González JM. Clinical Evaluation of Dental Implants with a Double Acid-Etched Surface Treatment: A Cohort Observational Study with Up to 10-Year Follow-Up. *Materials* (Basel). 2021 Oct 28;14(21):6483. doi: 10.3390/ma14216483. PMID: 34772010; PMCID: PMC8585230.

ATC (ASYMMETRICAL TRAPEZOIDAL COIL)

Kong L, Hu K, Li D, Song Y, Yang J, Wu Z, Liu B. Evaluation of the cylinder implant thread height and width: a 3-dimensional finite element analysis. *Int J Oral Maxillofac Implants*. 2008 Jan-Feb;23(1):65-74. PMID: 18416414.

McCullough JJ, Klokkevold PR. The effect of implant macro-thread design on implant stability in the early post-operative period: a randomized, controlled pilot study. *Clin Oral Implants Res*. 2017 Oct;28(10):1218-1226. doi: 10.1111/clr.12945. Epub 2016 Oct 3. PMID: 27699890.

Ao J, Li T, Liu Y, Ding Y, Wu G, Hu K, Kong L. Optimal design of thread height and width on an immediately loaded cylinder implant: a finite element analysis. *Comput Biol Med*. 2010 Aug;40(8):681-6. doi: 10.1016/j.compbiomed.2009.10.007. PMID: 20599193.

Abuhussein H, Pagni G, Rebaudi A, Wang HL. The effect of thread pattern upon implant osseointegration. *Clin Oral Implants Res*. 2010 Feb;21(2):129-36. doi: 10.1111/j.1600-0501.2009.01800.x. Epub 2009 Aug 25. PMID: 19709058.

Ryu HS, Namgung C, Lee JH, Lim YJ. The influence of thread geometry on implant osseointegration under immediate loading: a literature review. *J Adv Prosthodont*. 2014 Dec;6(6):547-54. doi: 10.4047/jap.2014.6.6.547. Epub 2014 Dec 17. PMID: 25551016; PMCID: PMC4279055.

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Certifications and Authorisations

EU-MDD
Certified

EU-MDR
Certified

Ente certificante | Certification Authority

Ente certificante | Certification Authority

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leadermedica.com

Registered office: Via Longhin, 11
35129 Padova Italy

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Leader Medica S.r.l. is responsible for both Certifications and Free Sale Certificates.

The other companies mentioned above are responsible for the world distribution of the devices made by Leader Medica S.r.l.

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Warning: All trademarks, trade names, and registered codes are the property of their respective owners and are mentioned exclusively to indicate compatibility with our devices.
Warning: Trademarks, trade names, and registered codes are the property of their respective owners and are mentioned exclusively to indicate compatibility with our devices.
LM Medical Technology (Shanghai) Co., Ltd
Sabettimed Srl
Tanto le Certificazioni che i Certificati di Libera Vendita sono in capo a Leader Medica S.r.l.