
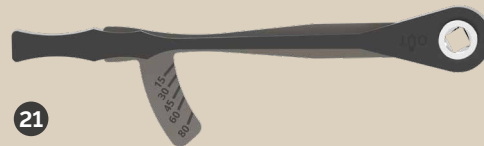
Kit Chirurgico MITO

MITO Surgical Kit

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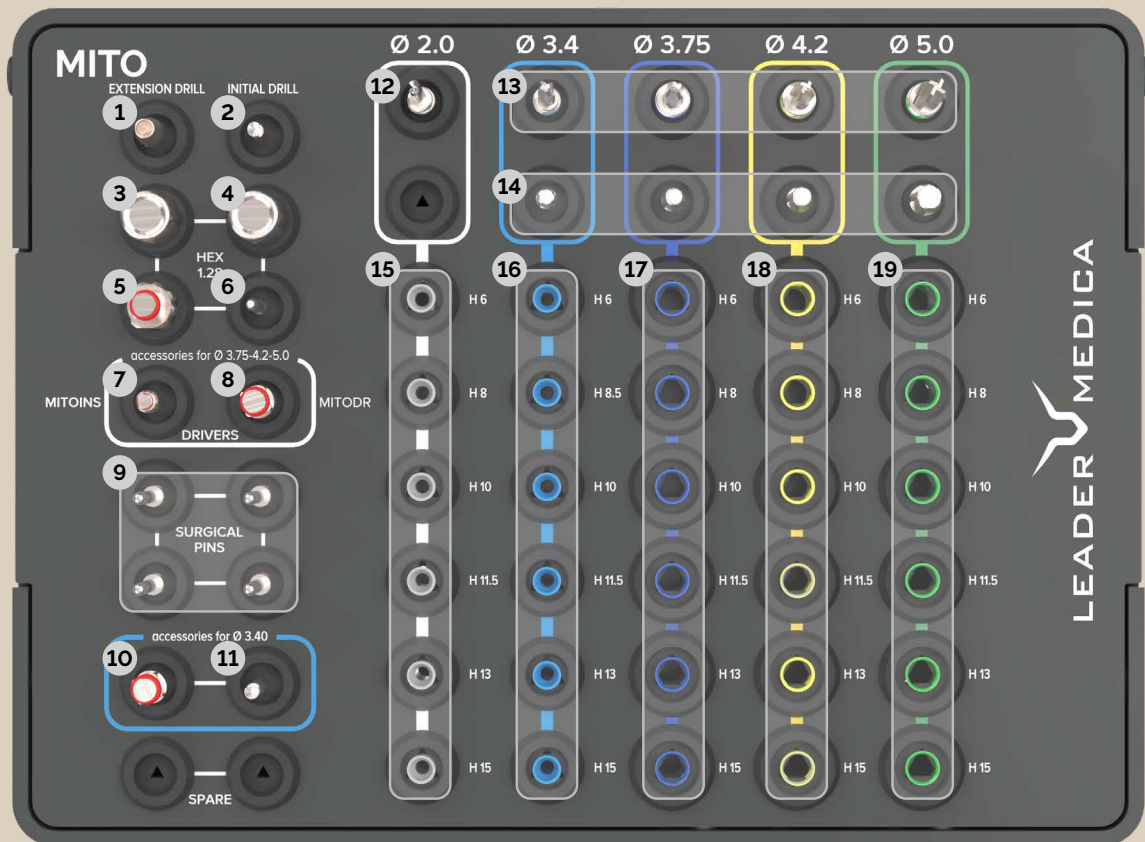For illustrative purposes, the surgical sequence for an implant with the following characteristics is shown below:

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
PLACEMENT OF THE SYSTEM AT CRESTAL LEVEL BONE LEVEL

For correct insertion, it is essential to check the crestal Bone Level positioning of the implant, in line with the bone crest (Fig. 1).

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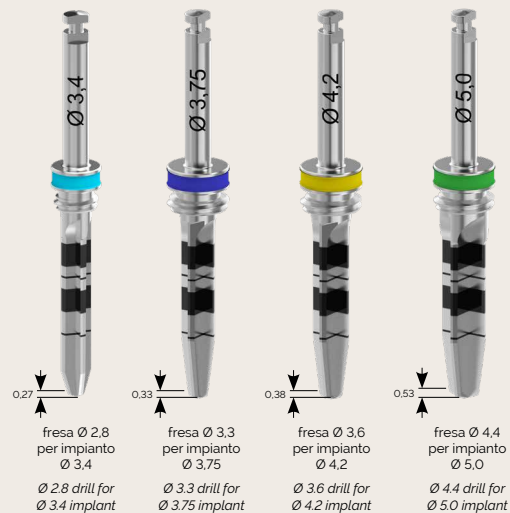IN ORDER TO DETERMINE THE EXACT MILLING DEPTH, IT IS RECOMMENDED TO CONSIDER THE LENGTH OF THE DRILL TIP, WHICH VARIES FOR EACH DRILL DIAMETER. (Fig. 2)

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The indications and images are just general instructions that should be adapted on a case-by-case basis according to the patient's needs. These recommendations are not intended to replace specialized training or the expertise of the clinician, who—also based on personal experience—can make different decisions.

ESSENTIAL PROCEDURES TO ENSURE THE SAFETY, EFFECTIVENESS, AND FUNCTIONALITY OF INSTRUMENTS

All surgical instruments must be thoroughly cleaned, disinfected, and sterilized before first use and before each subsequent reuse. It is essential to check the proper functioning of the instruments before each procedure and, if necessary, to immediately replace any instruments showing signs of wear or malfunction. Failure to follow these instructions may increase the risk of infection and intraoperative complications for the patient.

PREPARING THE AREA FOR IMPLANT PLACEMENT

Every procedure must be performed gradually and carefully, avoiding sudden movements that could compromise safety. For traditional surgery check the correct positioning of the drill stops, ensuring they are fully tightened. Choose the length in accordance with the implant insertion plan and the given warnings. Periodically check surgical micromotors used with drills and taps to prevent malfunctions. Regularly check the drills for wear and replace them if they lose effectiveness in order to ensure precision and safety during implant site preparation.

PROSTHESIS MAINTENANCE

To prevent potential mechanical and biological complications associated with prosthetic implants—which are widely documented in the literature—it is essential for the patient to maintain proper oral hygiene and undergo scheduled periodic checkups. These measures help prolong the functional life of the device. Periodic checks of the calibration of abutment and prosthetic screws are also recommended to ensure the stability and efficiency of the device. If a feeling of instability or malfunction occurs, the patient should be advised to seek prompt medical attention.

DISPOSING OF INSTRUMENTS AND IMPLANTS

Surgical instruments, consisting of small metal parts, can be disposed of as metal waste only if properly cleaned. Otherwise, they must be treated as biological waste. Similarly, dental implants removed from the patient's mouth must be treated as biological waste. For each disposal procedure, it is recommended to refer to the local regulations in force regarding the management of medical waste.